

ACTIVESHOP[®]

Instructions for use

Azzurro Lux cosmetic chair 3m



Importer:

ACTIVESHOP Sp. z o.o.
Ul. Graniczna 8b, building
DC2A 54-610 Wrocław,
Poland

www.activeshop.com.pl

Rev. A/07/2025



Scan for
different
languages



Attention!

1. The footrest is not designed to support the full weight of an adult.
2. Do not sit directly on the armrest and footrest.
3. Do not lift more weight than the chair can support.
4. For safety reasons, always plug the power cord into a grounded outlet.
5. Do not tie the power cord to the unit. Lifting or pulling the power cord may damage it. Make sure the power cord can be easily removed in an emergency.
6. Before moving the electric beauty chair, remove the power cord so that it does not get stuck between the parts of the chair.
7. If the power cord is damaged, immediately disconnect the cord, do not use the chair and contact an authorized service center. Use only the original power cord.
8. Do not let a person who is not performing the procedure control the control remote control.
9. Do not place any objects under the chair.
10. Make sure there is enough free space around the chair so that it can be folded and unfolded freely.
11. Make sure there are no sharp or harmful objects around the chair that could damage the device.
12. Do not modify the device yourself. If the chair does not work properly, stop using it and contact an authorized service center.
13. Do not throw any heavy objects on the chair or jump on it.

Mounting

1. Insert the headrest into the hole located at the top of the backrest (fig.1).
2. Position and connect the left and right armrest to the pin in the seat, and then rotate 90° counterclockwise (fig2.).

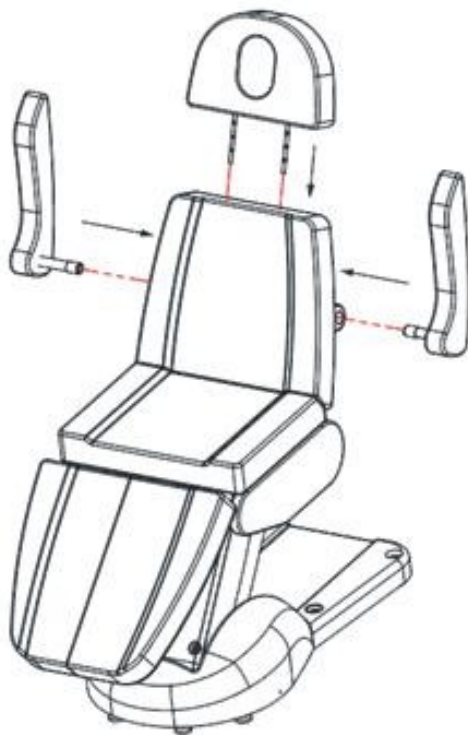


Fig.1



Fig.2


Remote control




1. Press the button  , to raise the seat.
2. Press the button  button to lower the seat.
3. Press the button  to increase the backrest angle.
4. Press the button  to decrease the backrest angle.
5. Press  to increase the angle of the footrest.
6. Press  to decrease the angle of the footrest.
7. Press  or  to turn on seat/back heating

Additional functions

- Pull the handle back and hold to adjust the angle of the footrest.

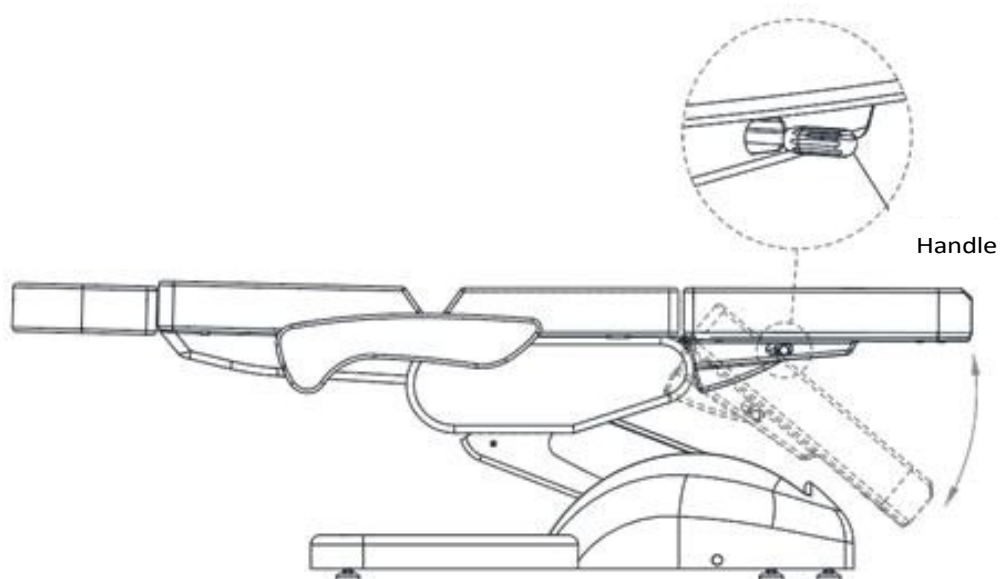

- The height of the headrest can be adjusted by pulling it up or pushing it down.

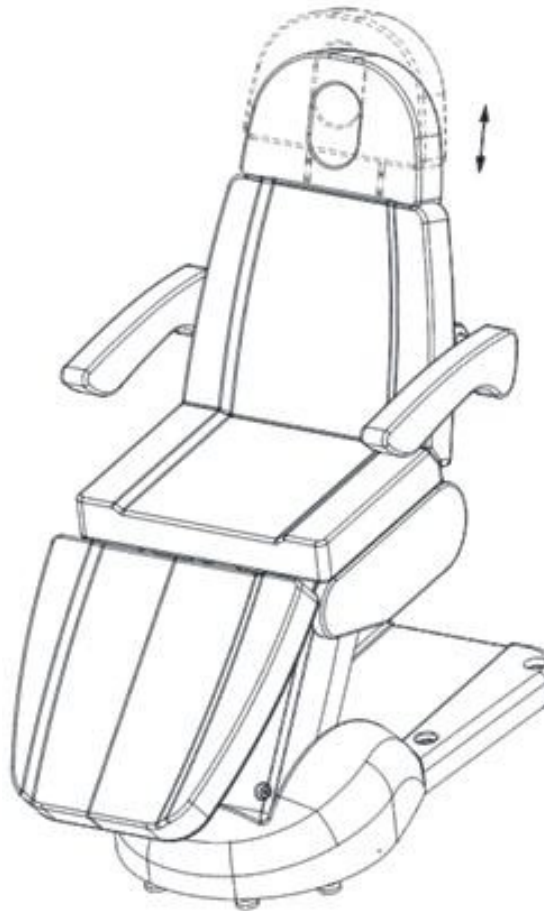

Maintenance measures:

Cleaning and disinfecting the device:

- Unplug the device before cleaning.
- To keep the surface of the chair in good condition, clean it regularly. Remember to clean and disinfect the chair after each use. It is recommended to perform a thorough disinfection once a month.

I Metal and plastic surfaces

- To clean metal and plastic surfaces and controls wipe them with a damp cloth. To clean hard-to-reach areas or corners, a small brush will work well. Do not use too much liquid to clean the device. Then disinfect the device with a suitable liquid.

II Stainless steel surfaces

- Cleaning with a mild detergent solution will work best. Ammonia and most solvents can be used to remove tough stains. Avoid chlorine-based solutions.

III Painted or chrome-plated metal surfaces

- Do not use water jets (showers, high-pressure water cannons) for cleaning.
- Do not clean in high temperature and high humidity, and do not expose the device to conditions that may lead to liquid flooding.
- Do not use solvents and acids for cleaning.

IV Chair surface

- Wipe the seat surface evenly with a special leather cleaner, then wipe it with a clean, soft cloth.

Unplug the device before maintenance.

- Repairs should be performed by experienced and manufacturer-authorized personnel.
- Repairs by unauthorized personnel may cause injury or damage to the device, for which the manufacturer is not responsible.
- Use only original spare parts approved by the manufacturer.
- Do not use the electric beauty chair or accessories if they are not working properly.
- All maintenance and repair work must be documented.

Attention!

- If you find any defects, stop using the equipment and contact an authorized service. Only experienced or authorized personnel can open or replace the actuator / control unit.
- If any part of the equipment is damaged, unplug the power cord and stop using it. Then notify the service department.

Troubleshooting:

Problem	Recognition	Solution
Buttons do not respond	<ul style="list-style-type: none">- Check that the power plug is connected- Check that all plugs connecting the transformer are properly connected	<ul style="list-style-type: none">- Disconnect the power plug and reconnect it- Disconnect all connecting plugs and reconnect them
There are noises and vibrations when lifting the chair	<ul style="list-style-type: none">- Loose screws	<ul style="list-style-type: none">- Tighten the screws
The indicator light is not on	<ul style="list-style-type: none">- The control remote control cable is not well connected or the remote control is damaged	<ul style="list-style-type: none">- Reconnect the cable, if it still does not work, replace the remote control

Specifications:

Dimensions: 2000 x 810 x 610-910

mm Maximum load: 140 kg Net

weight: 76 kg

Adjustable height range: 610-910 mm

Input voltage: AC 100-130 V / 220-240 V, 50/60 Hz Output

voltage: DC 24 V

SAFETY MEASURES

1. The device is a professional equipment for use in beauty parlors. The operation of the device only by qualified cosmetic personnel avoids undesirable consequences during treatment.

2. Do not disassemble the device or attempt to use it for purposes other than those indicated in this manual. All servicing and repairs must be performed only by an authorized service technician (ACTIV Tomasz Pacholczyk company).
3. Do not install or operate the device near liquids, wet places or with wet hands - so that no liquids get into the device.
4. Immediately unplug the device from the outlet and contact an authorized service provided by ACTIV Tomasz Pacholczyk in the following cases: a) liquid has gotten into the device, b) the device emits a strange smell, smoke or makes strange sounds, c) the power cord has been damaged, d) the device has been knocked over or its housing has been damaged.
5. Do not lay any objects on the cables. Do not install the device where someone may step on the power cord.
6. For safety purposes, unplug the power cord from the outlet each time you finish working with the device.
7. Do not insert anything into the opening of the device. If a foreign object gets into the hole, remove the cord from the outlet and contact the distributor or manufacturer.
8. Do not place the device on an unstable cart, shelf or vehicle. If dropped, the device may be damaged.
9. In order to guarantee the established useful life of the device, it is recommended to keep an interval of one minute before pressing the power button again (between switching on and off and vice versa).
10. The power supply must comply with the specifications on the nameplate. Otherwise, there is a risk of incorrect operation of the device its damage and even burning.

CAUTION:

1. In case of damage to the power cord or other malfunctions, contact an authorized service center - this will avoid danger.
2. This equipment may be used by children at least 8 years of age and by persons of diminished physical and mental capabilities and persons with inexperience and lack of familiarity with the equipment, if supervision or instruction on how to use the equipment in a safe manner is provided so that the associated risks are understood.
3. The equipment is not intended for play by children. Do not allow children to play with the equipment.
4. Do not leave the appliance unattended when children or persons without sufficient knowledge of the use of the appliance are in the vicinity.
5. This electrical appliance is not a toy. Improper use or contact with the electrical appliance may pose a health risk or lead to death, so use and store the appliance out of the reach of children.



DISPOSAL OF THE DEVICE

Dispose of used electrical appliances, fixtures and packaging to raw material reprocessing centers in accordance with environmental protection principles. Electrical appliances should not be disposed of in household waste containers. Follow the rules for discarding this type of equipment. In order to protect the environment in accordance with Directive 2012/19/EC concerning waste electrical and electronic equipment

Waste electrical and electronic equipment should be collected separately, and then handed over to collection points for waste electronic and electrical equipment.

At the end of their useful life, as well as old devices contain full-value recyclable components, so they should be handed over for reprocessing. Therefore, they should be handed over to appropriate collection points, as well as they can also be handed over to the ACTIV Tomasz Pacholczyk Company.

STORAGE AND TRANSPORT

The device should be transported at a temperature in the range of (-) 20 °C to (+) 50 °C only in its original packaging. Store the device in a well-ventilated and relative humidity of no more than 75%. Operating temperature from (+) 10 °C to (+) 40 °C. During storage and transportation processes, the device must not be thrown or dropped, exposed to mechanical shocks or atmospheric agents.

