

# ACTIVESHOP<sup>®</sup>

## USER MANUAL

Cosmetic chair hyd. Basic 210



### Importer:

ACTIVESHOP Sp. z o.o.  
Ul. Graniczna 8b, building  
DC2A 54-610 Wrocław,  
Poland

[www.activeshop.com.pl](http://www.activeshop.com.pl)

Rev. A/06/2024

Before using the cosmetic chair, read the enclosed instructions and perform the following steps:

1. Make sure all packaging items have been removed.
2. Make sure the chair is in working order.
3. Place the chair in the place where it will be used.

Attention.

Children and unauthorized persons should not use this device. Children must be supervised and must not play near the cosmetic chair.

### Dimensions of the device:



Length when fully unfolded: 184cm, Height min.


65cm,

Height max. 79cm

Maximum load of the chair: 150kg.



- The chair is manually adjustable with levers underneath. This way you can adjust the angle of the footrest and the backrest.
- The height of the chair is adjustable with the seat adjustment handle.
- In addition, the armrests are removable for easy access to the treatment areas.
- The cushion located on the headrest is also removable, making it easier to perform treatments when the person is lying on his or her stomach.
- The chair can be unfolded into a reclining position.
- The chair can be rotated 360°.

To adjust the position of the backrest/footrest, pull the lever located under the bottom of the chair and then tilt the backrest/footrest to the desired position.

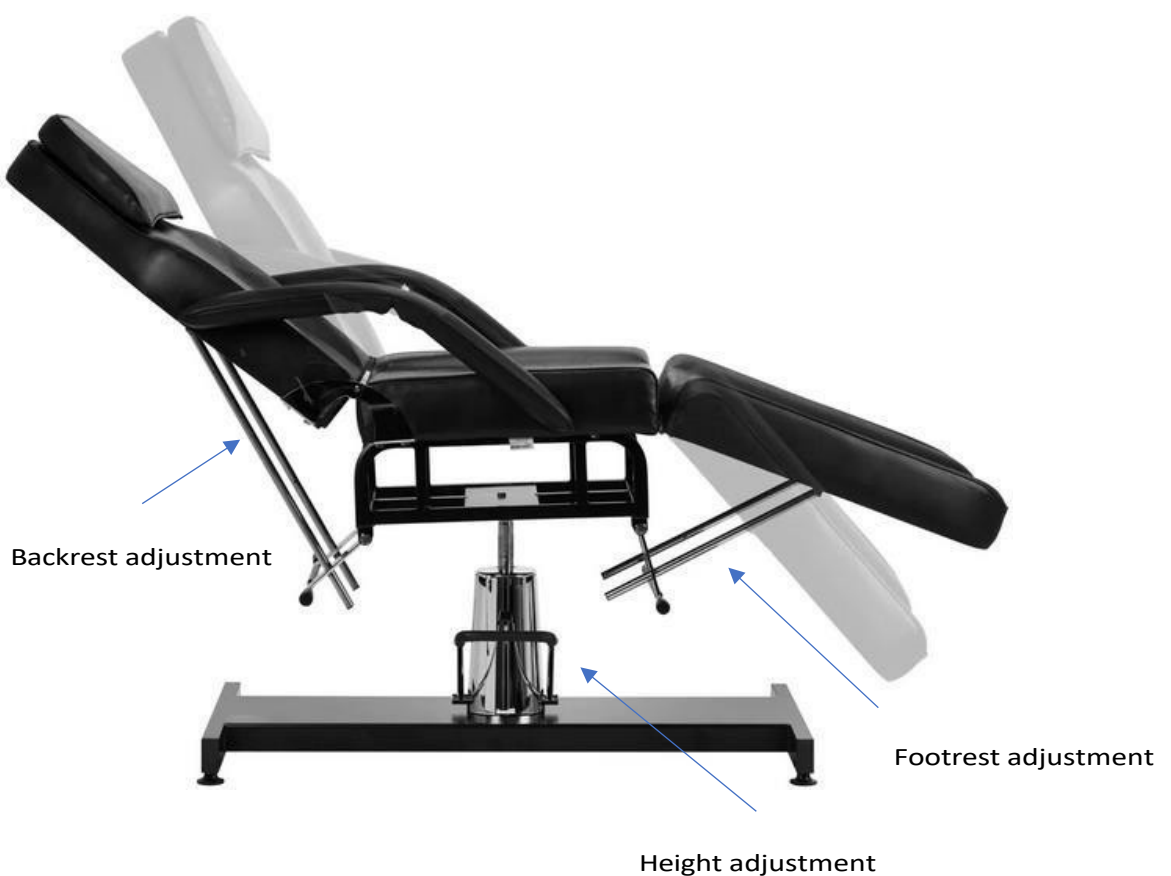

### **Assembly:**

1. Screw the seat to the frame.
2. Settle the seat on the actuator "press fit".
3. Tighten the actuator-seat locking screw.
4. Install the brackets in the backrest and seat.
5. Install the brackets in the footrests and seat.

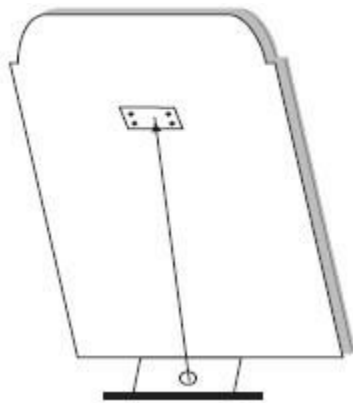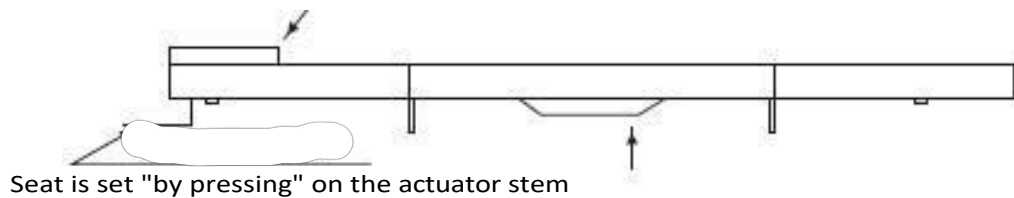

Pass the bracket through the hole, tighten with a screw.

In place on the backrest, we screw the bracket fastening it through the hole with a screw (the same for the footrests).

## **Cleaning and maintenance:**

- To keep the chair's surface in good condition, it should be cleaned regularly. Remember to clean and disinfect the chair after each use. It is recommended to perform a thorough disinfection once a month.

### **I Metal and plastic surfaces**

- To clean metal and plastic surfaces and controls, wipe them with a damp cloth. To clean hard-to-reach areas or corners, a small brush will work best.
- Do not use too much liquid to clean the device. Then disinfect the device with a suitable liquid.

### **II Stainless steel surfaces**

- Cleaning with a mild detergent solution will work best. Ammonia and most solvents can be used to remove tough stains. Chlorine-based solutions should not be avoided.

### **III Painted or chrome-plated metal surfaces**

- Do not use water jets (showers, high-pressure water cannons) for cleaning.
- Do not clean in high temperature and high humidity. Do not expose the device to conditions that may lead to liquid contamination.
- Do not use solvents and acids for cleaning.

#### **IV Chair surface**

- Wipe the seat surface evenly with a special leather cleaner, then wipe it with a clean, soft cloth.

#### **Maintenance**

##### **Unplug the device before maintenance.**

1. Repairs should be performed by experienced and manufacturer-authorized personnel. Repairs by unauthorized personnel may cause injury or damage to the device, for which the manufacturer is not responsible.
2. Use only original replacement parts approved by the manufacturer.
3. Do not use the electric beauty chair or accessories if they are not working properly.
4. All maintenance and repair work must be documented.

#### **Maintenance measures**

Check controls and their wiring once a week. Check the condition and function, of the following items at least once every six months:

1. Power cord and its attachment
2. Motor wiring
3. The control box and its wiring
4. Whether accessories are attached
5. Normal operation of the locking system

If you find any defects, stop using the equipment and contact an authorized service center. Only experienced or authorized personnel can open or replace the actuator / control unit.

If any part of the equipment is damaged, unplug the power cord and stop using it. Then notify the service center.

#### **SAFETY MEASURES**

1. The device is a professional equipment for use in beauty parlors. Operation of the device only by qualified cosmetic personnel avoids undesirable consequences during treatment.
2. Do not disassemble the device or attempt to use it for purposes other than those indicated in this manual. All servicing and repairs must be performed only by an authorized service technician (ACTIVESHOP L t d . company).
3. Do not install or operate the device near liquids, wet places or with wet hands - so that no liquids get into the device.
4. Immediately unplug the device from the outlet and contact the authorized service provided by ACTIVESHOP Ltd. in the following cases:  
(a) liquid has entered the device, (b) the device emits a strange odor, smoke or makes strange sounds,  
(c) the power cord has been damaged, (d) the device has been knocked over or its housing has been damaged.
5. Do not put any objects on the cables. Do not install the device where someone may step on the power cord.
6. For safety purposes, unplug the power cord from the outlet every time you finish working with the device.
7. Do not insert anything into the opening of the device. If a foreign object gets into the hole, remove the cord from the outlet and contact the distributor or manufacturer.
8. Do not place the device on an unstable cart, shelf or vehicle. If dropped, the device may be damaged.

9. In order to guarantee the established useful life of the device, it is recommended to maintain an interval of one minute before pressing the power button again (between switching on and off and vice versa).

10. The power supply must comply with the specifications on the nameplate. Otherwise, there is a risk of incorrect operation of the device its damage and even burning.

#### **CAUTION:**

1. In case of damage to the power cord or other malfunctions, contact an authorized service center - this will avoid danger.

2. This appliance can be used by children at least 8 years of age and by persons with reduced physical and mental capabilities and persons with inexperience and lack of familiarity with the equipment, if supervision or instruction on the use of the equipment in a safe manner is provided, so that the associated risks are understood.

3. The device is not intended for play by children. Do not allow children to play with the device.

4. Do not leave the appliance unattended when children or persons without sufficient knowledge in the use of the appliance are near it.

5. This electrical appliance is not a toy. Improper use or contact with the electrical appliance may pose a health risk or lead to death, so use and store the appliance out of the reach of children.

#### **DISPOSAL OF THE DEVICE**



Dispose of used electrical appliances, fixtures and packaging to raw material reprocessing centers in accordance with environmental protection principles. Electrical appliances should not be thrown into household waste containers. Follow the rules for discarding this type of equipment. In order to protect the environment in accordance with the directive

2012/19/EC concerning waste electrical and electronic equipment, waste electrical equipment should be collected separately, and then handed over to collection points for waste electronic and electrical equipment. At the end of their useful life, as well as old appliances contain full-value recyclable components, so they should be handed over for reprocessing. Therefore, they should be handed over to appropriate collection points, as well as they can also be handed over to ACTIVESHOP Ltd.

#### **STORAGE AND TRANSPORT**

The device should be transported at a temperature in the range of (-) 20 °C to (+) 50 °C only in its original packaging. Store the device in a well-ventilated and relative humidity of no more than 75%. Operating temperature from (+) 10 °C to (+) 40 °C. During storage and transportation processes, the device must not be thrown or dropped, exposed to mechanical shocks or atmospheric agents.