

ACTIVESHOP®

USER MANUAL

Azzurro Special 152 Beauty Chair



Importer:

ACTIVESHOP Sp. z o.o.
8b Graniczna St., Bldg. DC2A 54-610
Wrocław, Poland
www.activeshop.com.pl
Rev. A/02/2026

READ THE USER MANUAL BEFORE USING THE PRODUCT


This manual contains detailed information regarding the operation and maintenance of the device. Before using the device, please read the enclosed manual. Use the device as described. Keep this manual and ensure that all persons who may use the device are familiar with it.

Technical specifications:

Dimensions: 119 x 60 x 53 cm
Weight: 14 kg
Maximum weight capacity: 150 kg
Height adjustment: 55–75 cm Backrest
tilt adjustment: 0–5 cm Seat angle
adjustment: 10°

Construction:

1. Chair backrest and seat
2. Gas lift
3. Chair base
4. Casters

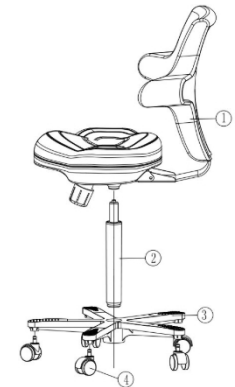

Assembly:

1. On some models, the chair backrest and seat are not pre-assembled. In this case, slide the backrest into the seat mounting and secure it with the 3 included screws.
2. Attach the casters to the base by pressing them into place.
3. Insert the gas lift with the thicker end into the base.
4. Place the seat onto the gas lift

Note: On some models, the gas spring is protected by a plastic cover (where it connects to the seat). Before installing the seat, make sure the plastic cover has been removed. Otherwise, it will not be possible to adjust the vertical position after assembly.

Operation:

Vertical adjustment:

To lower the seat, press the lever shown in the picture while sitting on the seat or otherwise applying weight to it.

To raise the seat, press the lever shown in the picture without applying weight to the seat at the same time.




Backrest tilt adjustment:

Pull the end of the lever outward to unlock the recline adjustment. Press the backrest and adjust it to the desired position.


Then push the end of the lever back in to lock the backrest in place.

Please note that due to the type of locking mechanism, to move the backrest closer to the seat, you must first push it in the opposite direction (as indicated by the arrow in attached illustration.)

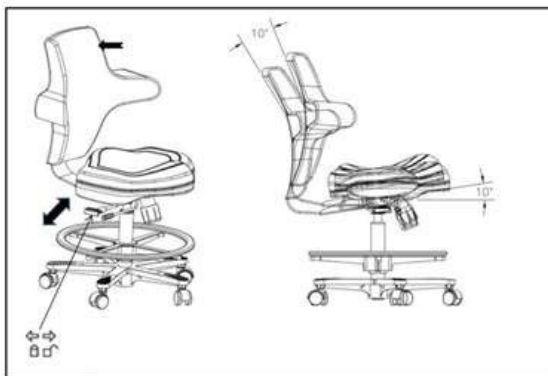

Horizontal seat adjustment:

To move the seat forward or backward, lift the marked in the picture, and then move the seat to the desired position. When you release the lever, the seat will remain in the selected position.

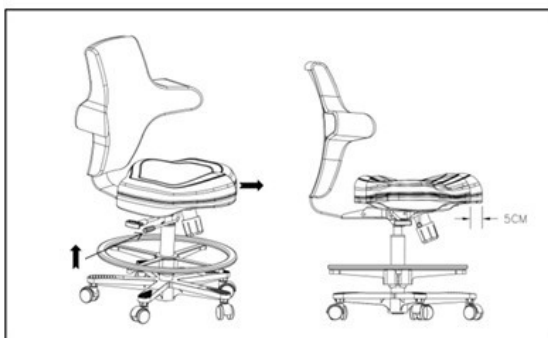

Recommendations:

1. Do not exceed the recommended weight limit for the chair; exceeding the recommended weight may cause permanent damage (max 150 kg).
2. Do not use products containing corrosive agents or dyes, or oily cleaning agents on the chair.
3. Inappropriate cleaning agents may cause damage to the material, corrosion, or other issues.
4. Clean the seat regularly to prevent dirt from becoming permanently embedded. If the upholstery gets dirty, it is safest to use lukewarm soapy water.

STORAGE AND TRANSPORT

The device should be transported at temperatures ranging from (–) 20 °C to (+) 50 °C exclusively in its original packaging. Store the device in a well-ventilated area with relative humidity not exceeding 75%. Operating temperature: from (+) 10 °C to (+) 40 °C. During storage and transport, do not throw or drop the device, or expose it to mechanical shocks or weather conditions.

SAFETY PRECAUTIONS

1. This device is professional equipment intended for use in beauty salons. To avoid adverse effects during treatment, the device must be operated exclusively by qualified beauty professionals.
2. Do not disassemble the device or attempt to use it for purposes other than those specified in this user manual. All maintenance and repairs must be performed exclusively by an authorized service technician (ACTIVESHOP Sp. z o.o.).
3. Do not install or operate the device near liquids, in wet areas, or with wet hands—ensuring that no liquids enter the device.
4. Immediately unplug the device from the outlet and contact the authorized service provided by ACTIVESHOP Sp. z o.o. in the following cases:
 - a) liquid has entered the device, b) the device emits a strange odor, smoke, or unusual noises, c) the power cord has been damaged, d) the device has been knocked over or its casing has been damaged.
5. Do not place any objects on the cables. Do not install the device in locations where someone might step on the power cord.
6. For safety reasons, unplug the power cord from the outlet after each use of the device.
7. Do not insert anything into the opening in the device. If a foreign object enters the opening, unplug the power cord and contact the distributor or manufacturer.
8. Do not place the device on an unstable cart, shelf, or vehicle. The device may be damaged if it falls.
9. To ensure the device's specified service life, it is recommended to wait one minute before pressing the power button again (between turning the device on and off, and vice versa).
10. The power supply must comply with the specifications listed on the nameplate. Otherwise, there is a risk of the device malfunctioning, becoming damaged, or even catching fire.

CAUTION:

1. In the event of damage to the power cord or other malfunctions, contact an authorized service center—this will help prevent hazards.
2. This device may be used by children aged 8 years and older and by persons with reduced physical or mental capabilities, as well as those with no experience or knowledge of the equipment, provided that supervision or instruction regarding the safe use of the equipment is provided so that the associated risks are understood.
3. This device is not intended for use by children. Do not allow children to play with the device.
4. Do not leave the device unattended when children or people who lack sufficient knowledge of how to use the device are nearby.
5. This electrical device is not a toy. Improper use or contact with the electrical device may pose a risk of injury or death; therefore, the device must be used and stored out of the reach of children.

DISPOSAL OF THE DEVICE



Used electrical devices, accessories, and packaging should be taken to recycling centers in accordance with environmental protection regulations. Electrical devices should not be disposed of in household waste bins. Follow the regulations regarding disposal of this type of equipment. To protect the environment in accordance with Directive 2012/19/EU on waste electrical and electronic equipment, used electrical devices must be collected separately and then taken to collection points for waste electrical and electronic equipment. At the end of their useful life, as well as older devices, contain valuable components suitable for recycling, so they should be sent for recycling. Therefore, they should be taken to appropriate collection points, and they can also be sent to ACTIVESHOP Sp. z o.o.