

ACTIVESHOP[®]


INSTRUCTIONS FOR USE

Azzurro 891 electric cosmetic chair



Importer:

ACTIVESHOP Sp. z o.o.
Ul. Graniczna 8b, building
DC2A 54-610 Wrocław,
Poland
www.activeshop.com.pl
Rev. A/06/2024

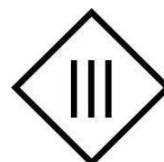


Before using the cosmetic chair, read the instructions and perform the following steps:

1. Make sure all packaging materials have been removed.
2. Make sure the chair is in working order.
3. Place the cosmetic chair in the area where it will be used.
4. Plug the power plug into a suitable outlet.

Tip! Children and inexperienced people must not use this device. Children must be supervised and must not play near the cosmetic chair.

Dimensions of the device:

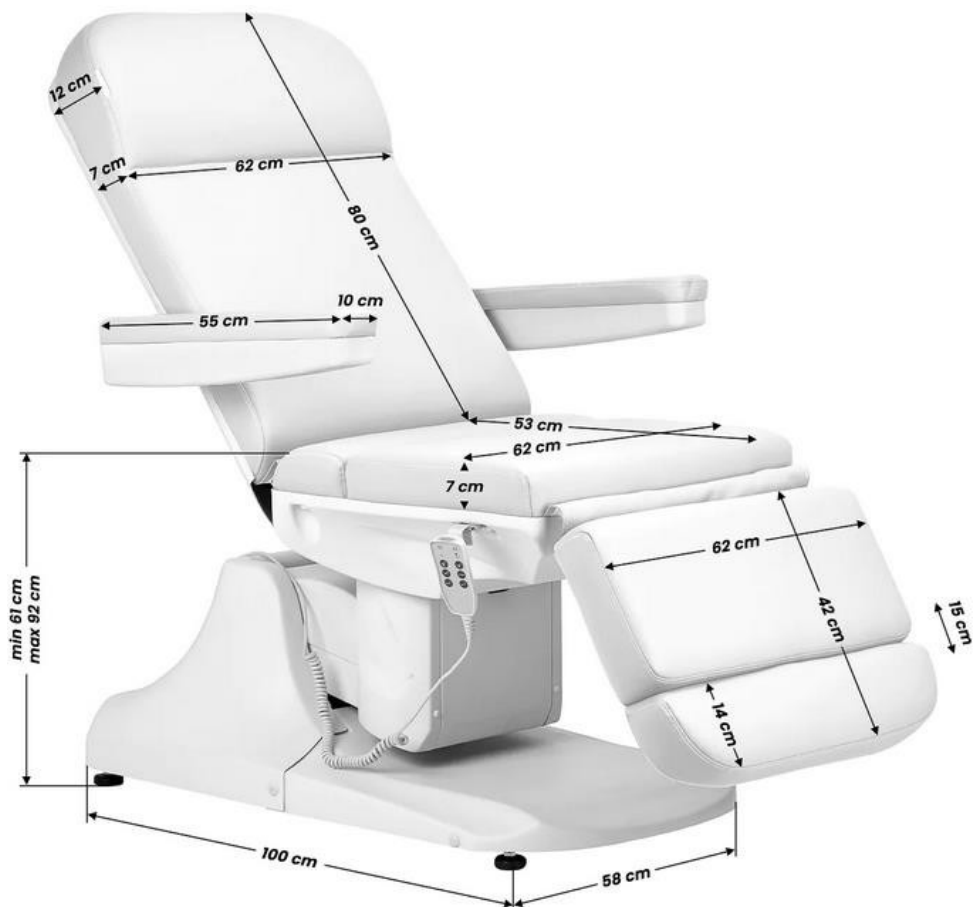

The device is controlled by a remote control.

Functions of the remote control:



1. Seat height adjustment from 61-92cm (press to raise or lower the seat).
2. Backrest angle adjustment
3. Adjustment of the footrest (press to change the position of the footrest)

Footrest adjustment:

The footrest has an adjustment option and can be extended from 0-15cm. To extend the footrest, turn the knob on the bottom and extend the footrest to the desired length.



Cleaning and disinfection of the device

Disconnect the device from the power supply before cleaning.

To keep the chair surface in good condition, clean it regularly. Remember to clean and disinfect the chair after each use. It is recommended to perform a thorough disinfection once a month.

I Metal and plastic surfaces

To clean metal and plastic surfaces and controls wipe them with a damp cloth. To clean hard-to-reach areas or corners, a small brush will work well. Do not use too much liquid to clean the device. Then disinfect the device with a suitable liquid.

II Stainless steel surfaces

Best will work will be cleaning with a mild solution detergent. Ammonia i most solvents can be used to remove tough stains. Avoid chlorine-based solutions.

III Painted or chrome-plated metal surfaces

- Do not use water jets (showers, high-pressure water cannons) for cleaning.
- Do not clean in high temperature and high humidity, and do not expose the device to conditions that may lead to liquid contamination.
- Do not use solvents and acids for cleaning.

IV Seat surface

- Wipe the surface of the seat evenly with a special leather cleaner, then wipe it with a clean, soft cloth.

Maintenance

Unplug the device before maintenance.

1. Repairs should be performed by experienced and manufacturer-authorized personnel. Repairs by unauthorized personnel may cause injury or damage to the device, for which the manufacturer is not responsible.
2. Use only original spare parts approved by the manufacturer.
3. Do not use the electric beauty chair or accessories if they are not working properly.
4. All maintenance and repair work must be documented.

Maintenance measures

Check the controls and their wiring once a week. Check the condition and function, of the following items at least once every six months:

1. Power cord and its attachment
2. Motor wiring
3. The control box and its wiring
4. Whether accessories are attached
5. Normal operation of the locking system

If you find any defects, stop using the equipment and contact an authorized service center. Only experienced or authorized personnel can open or replace the actuator / control unit.

If any part of the equipment is damaged, unplug the power cord and stop using it. Then notify the service center.

Troubleshooting

Problem	Cause	Solution method
When you press a button on the remote control, the remote control does not respond	1. the function of a particular button has already reached its maximum.	1. try to press another button
The device makes strange sounds when adjusting height	Loose screw	Tighten all screws