

ACTIVESHOP[®]

USER MANUAL

**Elegante 6027 Magnifying
Lamp, 60 SMD LEDs, 5D,
for tabletops**



Importer:
ACTIVESHOP SP. Z O.O.
8b Graniczna St., Bldg. DC2A
54-610 Wrocław, Poland
www.activeshop.com.pl
Rev. A/10/2025



Scan for
different
languages



Thank you for choosing our product.
Please read the user manual before first use and keep it for future reference.

The Elegante 6027 Magnifying Lamp with tabletop mount features high-quality construction and a very stable arm. The two-piece arm offers a wide range of tilt adjustment, allowing you to perfectly adjust the light direction to your needs. Thanks to SMD LEDs, the lamp provides bright light while consuming minimal power. The glass lens under the flap allows you to view the skin under magnification while providing excellent illumination. Locking screws ensure a secure and stable position of the arm in the selected position.


Technical Specifications:	
Number of LEDs	60 LEDs
Light output	820 lumens
Color temperature	6,500 K
Voltage	220–240 V, 50 Hz
Power	10W

Installation Instructions:

- 1. Screw the mounting bracket to the tabletop
- 2. Insert the lamp post into the mounting bracket
- 3. Plug the cord into a power outlet
- 4. Adjust the lamp's angle and lock it in place using the locking screw
- 5. The power button is located on the top of the lamp
- 6. Lift the flap to use the magnifying glass

CAUTION - Do not expose unprotected eyes to direct LED light, as this may cause discomfort or damage to your eyesight.

Product Dimensions

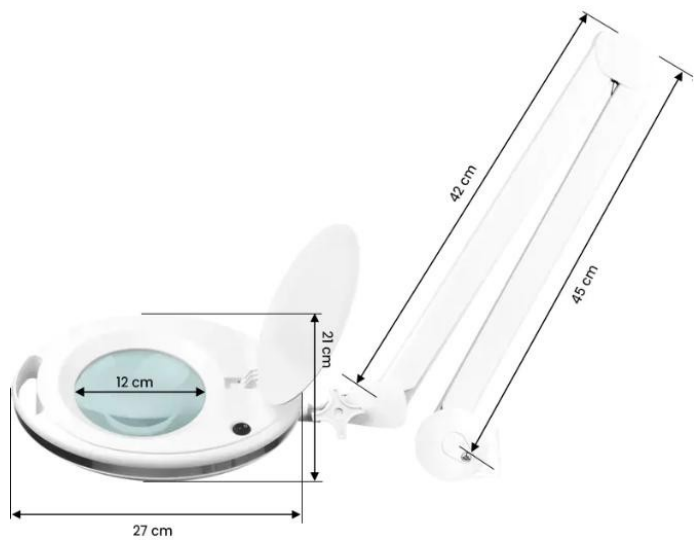

Maintenance

To keep the product in good working condition, ensure a long service life, and maintain proper operation, we recommend regular maintenance.

- 1. Before cleaning, turn off the device, unplug it from the outlet, and make sure the lamp has cooled down.
- 2. Wipe the device with a soft cloth.

3. Do not use harsh detergents that could damage or scratch the surface.

Recommendations

1. Make sure you place the device on a flat, stable surface to reduce the risk of it falling and becoming damaged.
2. Never leave the device plugged in and unattended.
3. Use the device as intended, as described in the manual.
4. Never use this device if the cord or plug is damaged.
5. Do not use the device if it is not working properly or if it has been dropped, damaged, or submerged in water.
6. Do not pull directly on the cord to unplug the device from the outlet.
7. This product is for indoor use only. Do not expose it to direct weather conditions.
8. Clean the device periodically to prevent dirt buildup. It is best to clean the device with a damp cloth.
9. Do not wrap the power cord around the device.

Safety Precautions

1. This device is professional equipment intended for use in beauty salons. To avoid adverse effects during treatment, the device must be operated exclusively by qualified beauty professionals.
2. Do not disassemble the device or attempt to use it for purposes other than those specified in this user manual. All maintenance and repairs must be performed exclusively by an authorized service technician (ACTIVESHOP Sp. z o.o.).
3. Do not install or operate the device near liquids, in wet areas, or with wet hands. Ensure that no liquids enter the device.
4. Immediately unplug the device from the outlet and contact the authorized service center provided by ACTIVESHOP Sp. z o.o. in the following cases: a) liquid has entered the device, b) the device emits a strange odor, smoke, or makes unusual noises, c) the power cord has been damaged, d) the device has been knocked over or its housing has been damaged.
5. Do not place any objects on the cables. Do not install the device in locations where someone might step on the power cord.
6. For safety reasons, unplug the power cord from the outlet after each use of the device.
7. Do not place any foreign objects inside the device. If a foreign object gets inside, unplug the power cord from the outlet and contact the distributor or manufacturer.
8. Do not place the device on an unstable cart, shelf, or vehicle. If it falls, the device may be damaged.
9. To ensure the device's specified service life, it is recommended to wait one minute before pressing the power button again (between turning it on and off, and vice versa).
10. The power supply must comply with the specifications listed on the nameplate. Otherwise, there is a risk of the device malfunctioning, being damaged, or even catching fire.

CAUTION

1. In the event of damage to the power cord or other malfunctions, contact an authorized service center—this will help prevent hazards.
2. This device may be used by children aged 8 and older, by persons with reduced physical or mental capabilities, and by persons unfamiliar with how this equipment, provided they are given appropriate supervision or instruction on how to use it safely and understand the associated risks.
3. Do not leave the device unattended when children or people without sufficient knowledge

regarding the use of the device are in its vicinity.

4. The device is not intended for use as a toy by children. Do not allow children to play with the device.
5. This electrical device is not a toy. Improper use of or contact with the device
Electrical hazards may pose a risk to health or lead to death; therefore, the device should be used and stored out of the reach of children.

Disposal of the Device



To protect the environment, in accordance with Directive 2012/19/EC on waste electrical and electronic equipment, used electrical devices must be collected separately and taken to collection points for waste electrical and electrical equipment. Follow the guidelines for disposing of this type of equipment and do not place it in household waste bins.

At the end of their useful life, these devices still contain fully functional components that can be recycled; therefore, they should be sent for recycling to appropriate collection points or to ACTIVESHOP Sp. z o.o.

Storage and Transport

The device should be transported at temperatures ranging from (–)20°C to (+)40°C exclusively in its original packaging. It should be stored in a well-ventilated room with relative humidity not exceeding 90%. Operating temperature range: (+)10°C to (+)40°C. During storage and transport, do not throw or drop the device, or expose it to mechanical shocks or weather conditions.