

ACTIVESHOP[®]

USER MANUAL

Podiatrist's Assistant 984

HZ-984 gray / white



Importer:
ACTIVESHOP SP. Z O.O.
8b Graniczna St., Bldg. DC2A
54-610 Wrocław, Poland
www.activeshop.com.pl
Rev. A/06/2026



Scan for
different
languages



Thank you for choosing our product.

Please read the user manual before first use and keep it for future reference.


The 984 Podiatry Assistant is a simple and functional piece of furniture designed for podiatry clinics, beauty salons, and other treatment areas where good workplace organization is essential. It facilitates the storage of tools, cosmetics, and equipment, bringing order and aesthetics to the interior.

Technical Specifications	
Ultraviolet lamp power	6 W
Voltage	230 V ~ / 50 Hz
Maximum load capacity	30 kg

Integrated UV sterilizer in the drawer

The drawer with a built-in UV lamp ensures hygienic protection of instruments between procedures. The mirrored interior finish promotes even distribution of UV radiation, increasing the effectiveness of sterilization.

- The system is activated only when the drawer is fully closed,
- Pulling out the drawer immediately turns off the UV lamp.

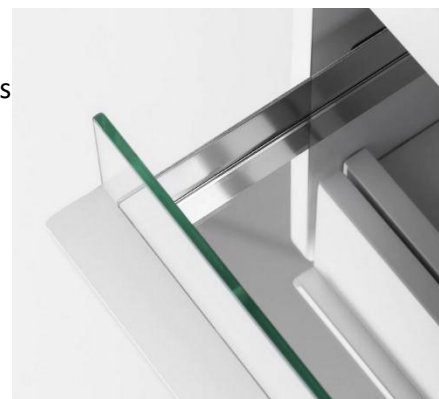

Integrated power outlets

The unit features a rear panel with three power outlets, two USB ports, and one USB-C port, allowing for convenient connection of various devices.



Ergonomic design

- Cable grommets located in the side panels near the top drawer.
- An integrated magnifying lamp holder allows for stable mounting of the lamp arm close to the treatment area.
- Practical casters facilitate easy movement of the cabinet within the office.

Product maintenance and care

- Light surface dirt and dust should be removed regularly with a soft, dry cloth.
- For greasy or stubborn stains, first wipe the cabinet with a cloth soaked in a mild detergent solution, then wipe the entire surface dry.
- Do not use strong solvents such as benzene, gasoline, or paint thinners, as they can irreversibly damage the paint finish.
- Do not expose the product to direct sunlight or store it in rooms with excessive humidity.

Product Dimensions



Recommendations

1. Make sure you install the device on a flat, stable surface to reduce the risk of it falling and being damaged.
2. Never leave the device plugged in and unattended.
3. Use the device as intended, as described in the manual.
4. Never use this device if the cord or plug is damaged.
5. Do not use the device if it is not working properly or if it has been dropped, damaged, or has come into contact with water.
6. Do not pull directly on the cord to remove the plug from the outlet.
7. This product is for indoor use only. Do not expose it to direct weather conditions.
8. Clean the device periodically to prevent dirt buildup. It is best to clean the device with a damp cloth.
9. Do not wrap the power cord around the device.

Troubleshooting

If the cabinet malfunctions, perform the following troubleshooting steps before contacting customer service:

Symptom	Probable cause	Remedy
Unusual noises are coming from inside the cabinet, and the structure is vibrating noticeably.	Loosening of mounting screws during operation.	Carefully check and tighten all screws connecting the cabinet's structural components.

Safety precautions

1. This device is professional equipment intended for use in beauty salons. Operation of the device exclusively by qualified beauty professionals helps prevent adverse effects during treatment.
2. Do not disassemble the device or attempt to use it for purposes other than those specified in this user manual. All service and repair work must be performed exclusively by an authorized service technician (ACTIVESHOP Sp. z o.o.).
3. Do not install or operate the device near liquids, in wet areas, or with wet hands. Ensure that no liquids come into contact with the device.
4. Immediately unplug the device and contact the authorized service provided by ACTIVESHOP Sp. z o.o. in the following cases: a) liquid has entered the device, b) the device emits a strange odor, smoke, or makes unusual noises, c) the power cord has been damaged, d) the device has been knocked over or its housing has been damaged.
5. Do not place any objects on the cables. Do not install the device in locations where someone might step on the power cord.
6. For safety reasons, unplug the power cord from the outlet after each use of the device.
7. Do not place foreign objects inside the device. If a foreign object gets inside, unplug the power cord from the outlet and contact the distributor or manufacturer.
8. Do not place the device on an unstable cart, shelf, or vehicle. The device may be damaged if it falls.
9. To ensure the device's specified service life, it is recommended to wait one minute before pressing the power button again (between turning it on and off and vice versa).
10. The power supply must comply with the specifications listed on the nameplate. Otherwise, there is a risk of the device malfunctioning, being damaged, or even catching fire.

CAUTION

1. In the event of damage to the power cord or other malfunctions, contact an authorized service center—this will help prevent hazards.
2. This device may be used by children aged 8 years and older, by persons with reduced physical or mental capabilities, as well as by persons unfamiliar with the operation of this equipment, provided they are given appropriate supervision or instruction regarding the correct and safe use of the device and understand the associated risks.
3. Do not leave the device unattended when children or persons without sufficient knowledge of how to use the device are in the vicinity.
4. The device is not intended for use by children. Do not allow children to play with the device.

5. This electrical device is not a toy. Improper use of or contact with the may pose a risk of injury or death; therefore, the device must be used and stored out of the reach of children.

Disposal of the device



To protect the environment, in accordance with Directive 2012/19/EC on waste electrical and electronic equipment, used electrical devices must be collected separately and taken to collection points for waste electronic and electrical equipment. Follow the rules for disposing of this type of equipment and do not dispose of it in household waste bins.

At the end of their useful life, these devices still contain valuable components suitable for recycling; therefore, they should be sent for recycling to appropriate collection points or to ACTIVESHOP Sp. z o.o.

Storage and Transport

The device should be transported at temperatures ranging from (–)20°C to (+)40°C exclusively in its original packaging. Store it in a well-ventilated room with relative humidity not exceeding 90%. Operating temperature range: (+)10°C to (+)40°C. During storage and transport, the device must not be thrown or dropped, nor exposed to mechanical shocks or weather conditions.