

ACTIVESHOP[®]

USER MANUAL

Cosmetic chair hyd. Sillon 210



Importer:

ACTIVESHOP Sp. z o.o.

ul. Graniczna 8b, bud. DC2A

54-610 Wrocław, Poland

www.activeshop.com.pl

Rev. A/05/2025

Scan for
different
languages




Dimensions:





Assembly:

A 4pcs	B M6 4pcs	C M6X35 4pcs	D 400 2pcs	E M5X8 8pcs	F 4pcs	G 250 2pcs	H M12X80 4pcs
I Ø12 4pcs	J M12 4pcs	K 500 2pcs	L M8X25 6pcs	M Ø8 6pcs	N Ø8 6pcs	O M10X25 1pcs	P Ø10 1pcs
Q Ø10 1pcs	R 250 2pcs	S 1pcs	T M6 4pcs	U M6 4pcs	V M6 4pcs	① 	
④ 						② 	
						③ 	
						CH-210	


1. Slide the pads into the base of the chair.
2. Mount the height adjustment lever on the actuator using the R element.
3. Settle the actuator on the seat base, then tighten the actuator using the V, U and T elements.
4. Place the metal seat frame on the base of the chair, and then fasten it, using elements L, M and N.

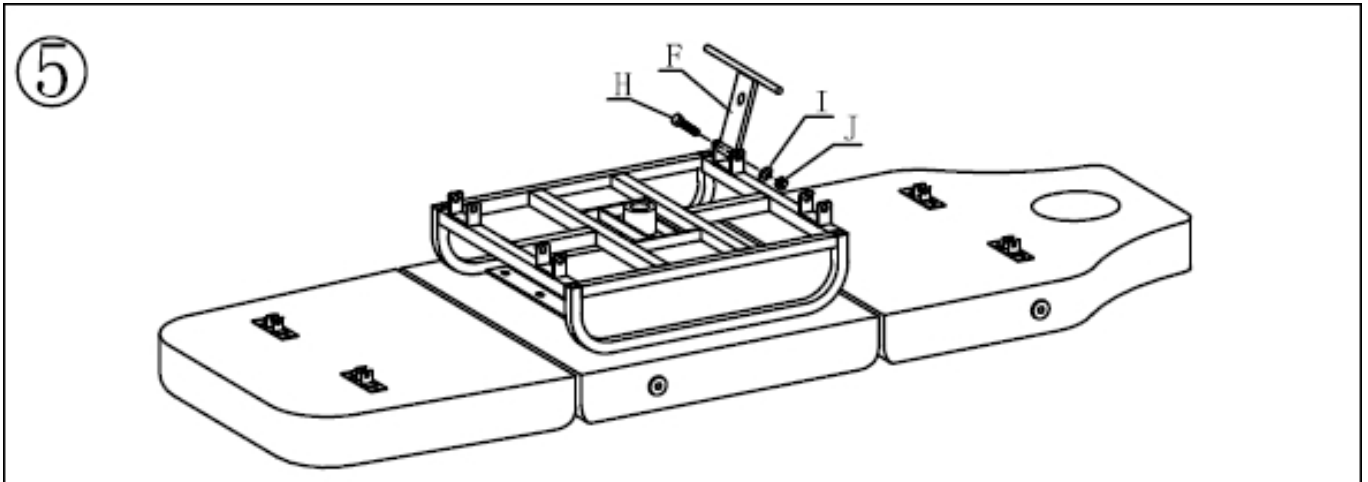


5. Install element F by tightening it with elements I, J and H.

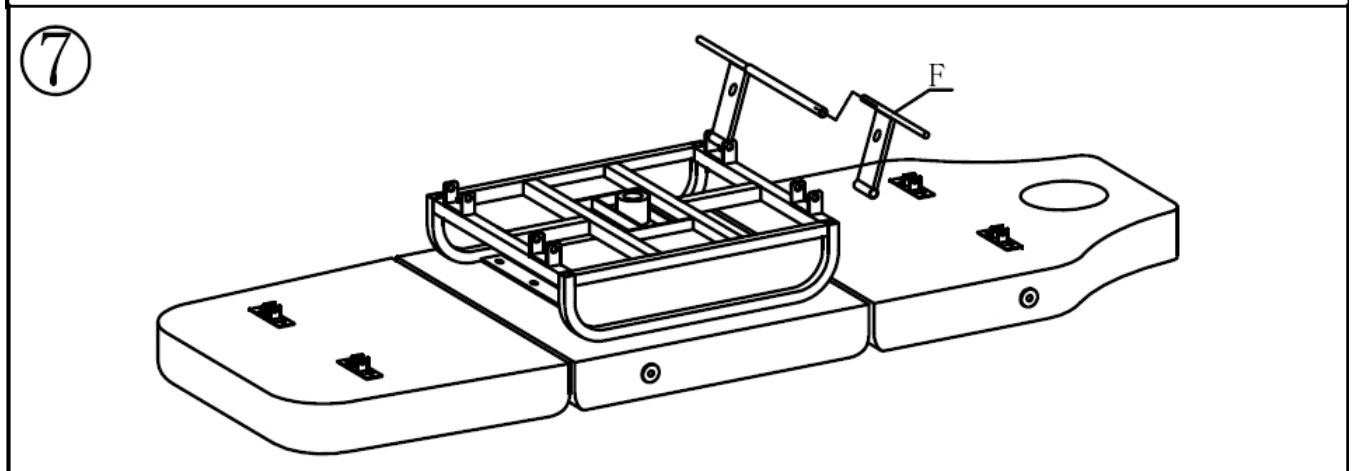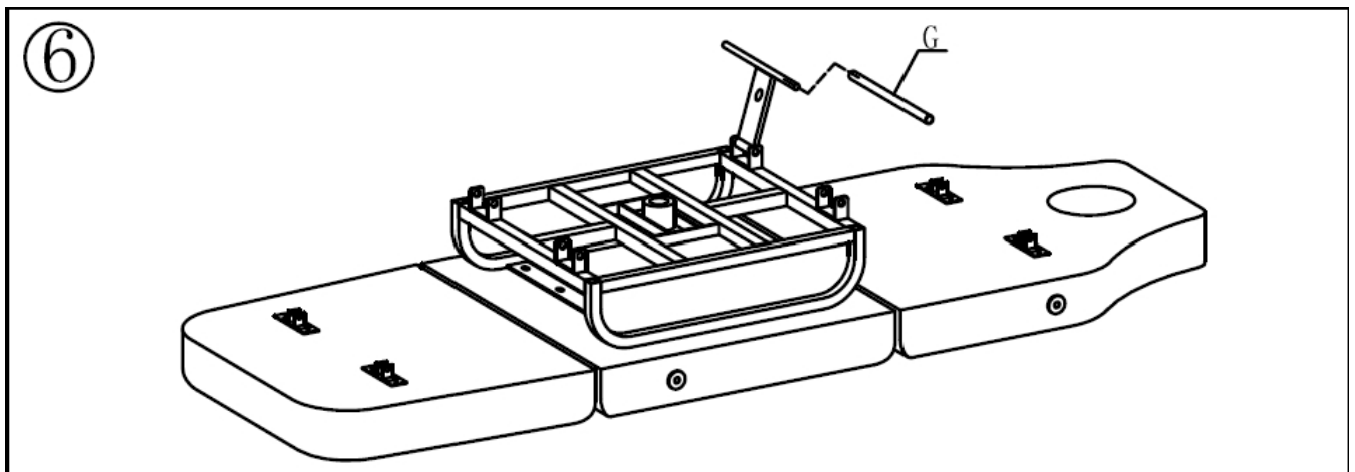


- 6,7. Use element G to connect the two F elements together.

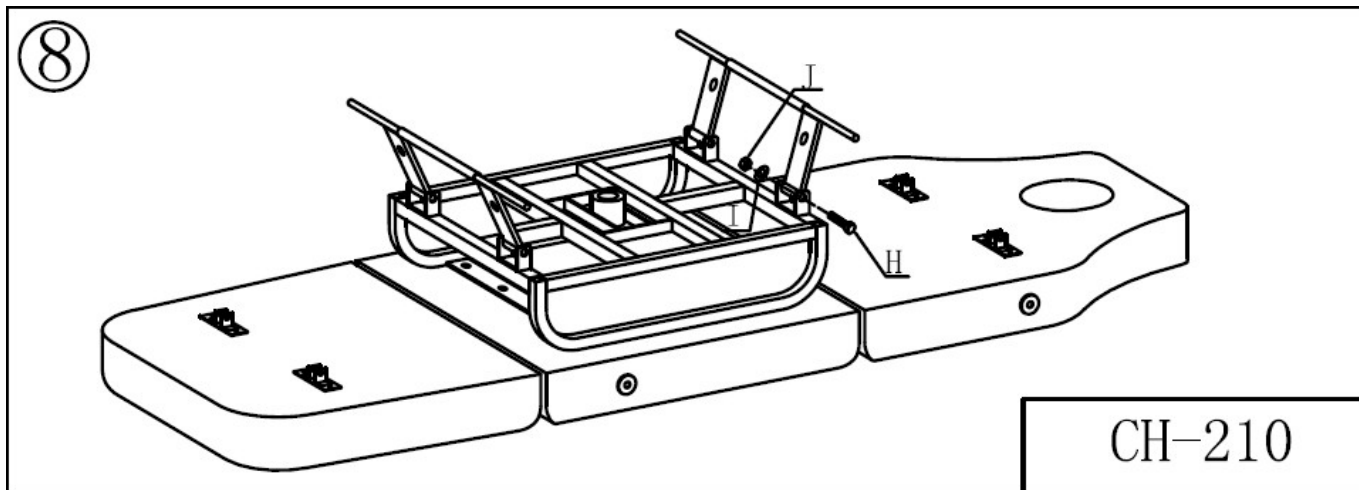


8. Tighten, using elements J and H, and then repeat on the other side.

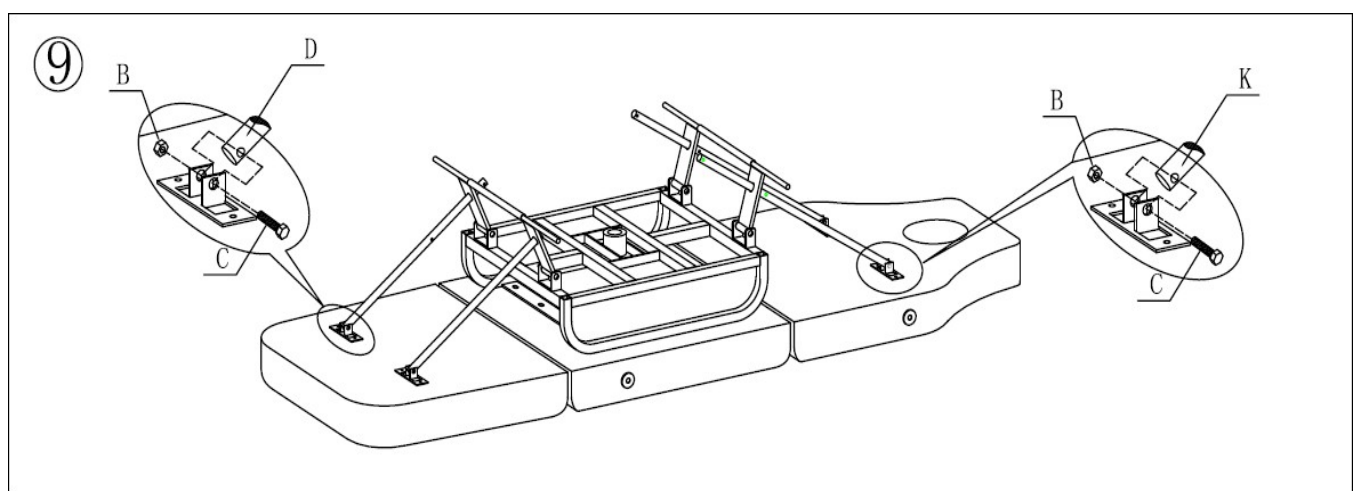

9. Insert elements D and K, and then screw them in place using B and C.

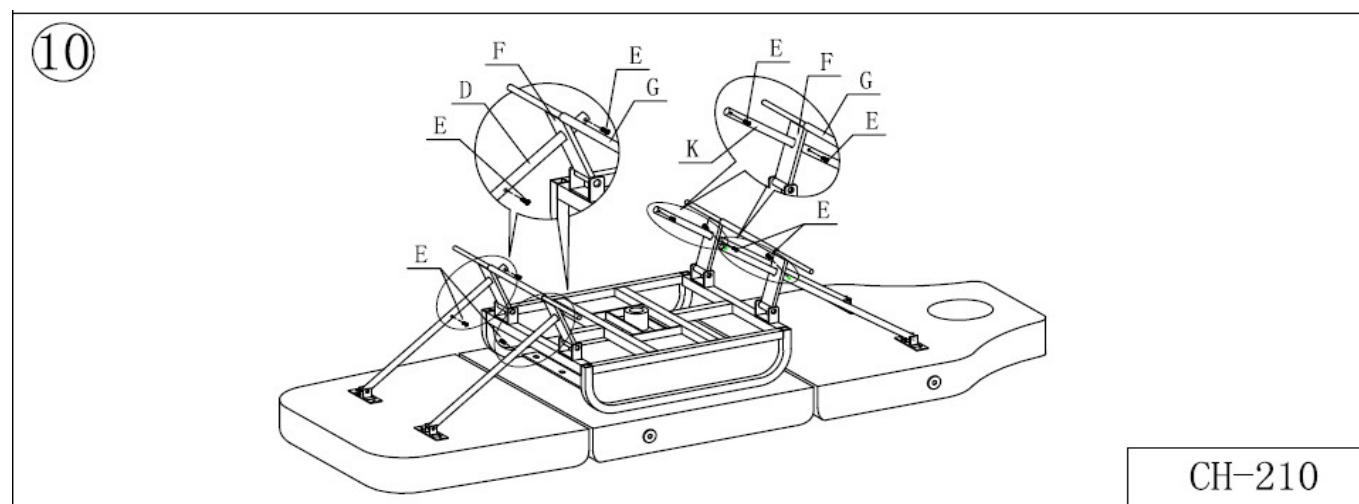

10. Secure the supports, using elements E.

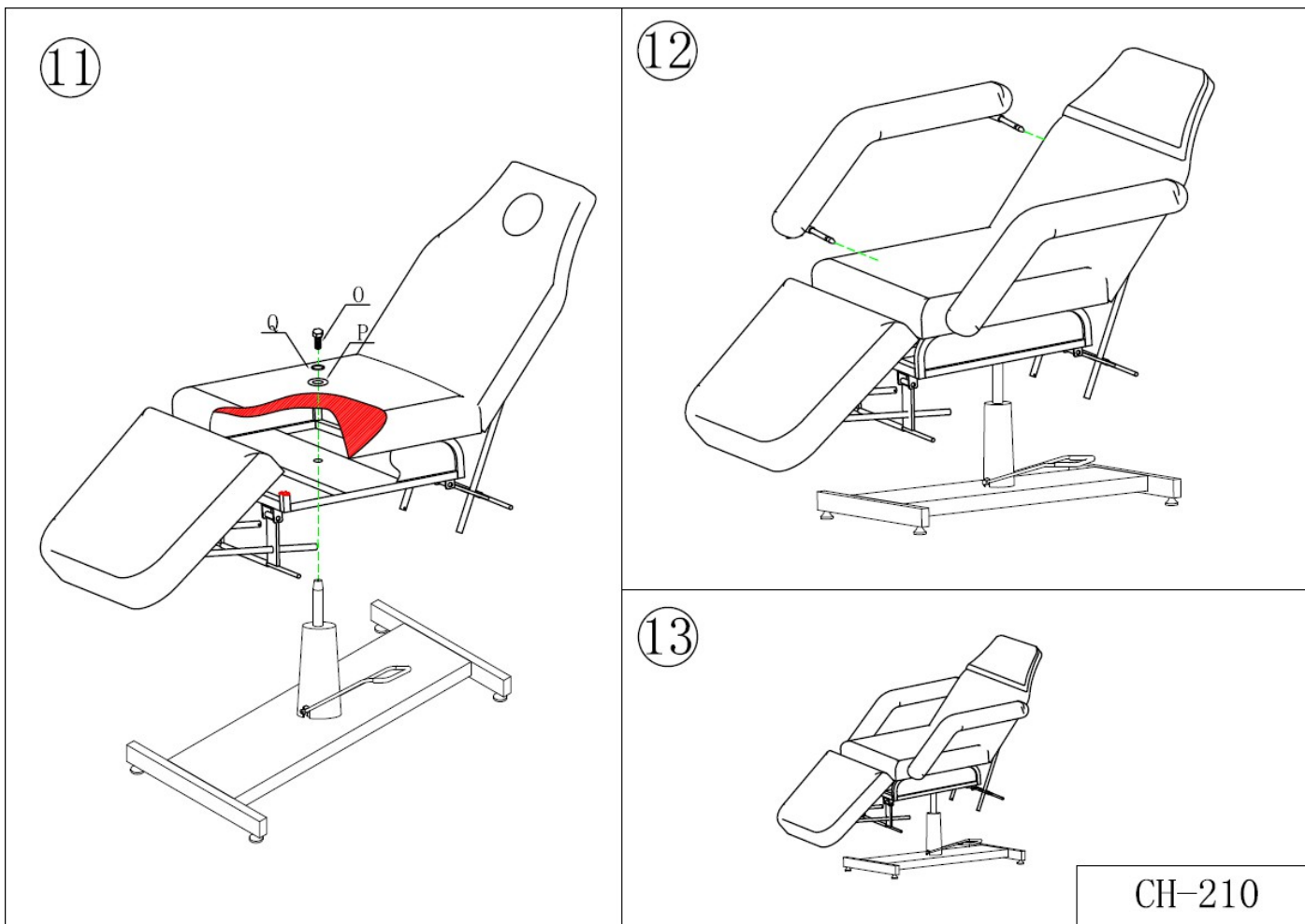

11. Set the seat on the actuator stem, then tighten it, using elements Q, P and O.

12. Slide the armrests in as shown in the illustration.

Seat functions:

- Lifting the seat: raise the seat by pressing down on the lever
- Lowering the seat: push down until it drops completely
- Locking the seat rotation: lift the lever as far up as possible

Caution.

Do not sit down abruptly or throw loads on the seat. It is forbidden to move or carry the chair while people or other loads are on it. Do not stand on the chair with your feet or knees, swing, sit on the armrests, and do not sit more than one person.

Cleaning and disinfecting the device:

Unplug the device before cleaning.

To keep the chair's surface in good condition, clean it regularly. Remember to clean and disinfect the chair after each use. It is recommended to perform a thorough disinfection once a month.

I Metal and plastic surfaces

To clean metal and plastic surfaces and controls wipe them with a damp cloth. To clean hard-to-reach areas or corners, a small brush will work well. Do not use too much liquid to clean the device. Then disinfect the device with a suitable liquid.

II Stainless steel surfaces

Cleaning with a mild detergent solution will work best. Ammonia and most solvents can be used to remove tough stains. Avoid chlorine-based solutions.

III Painted or chrome-plated metal surfaces

- Do not use water jets (shower, high-pressure water cannons) for cleaning.
- Do not clean in high temperature and high humidity, and do not expose the device to conditions that may lead to liquid contamination.
- Do not use solvents and acids for cleaning.

IV Chair surface

- Wipe the seat surface evenly with a special leather cleaner, then wipe it with a clean, soft cloth.

Maintenance:

1. **Daily cleaning**
Remove dust with a soft, dry cloth.
2. **Removing stubborn dirt**
As you know, it is difficult to remove stubborn dirt with an ordinary dry cloth. In this case, first use a soft cloth with diluted neutral detergent, and then wipe with a dry cloth.
3. **The use of benzene, thinners, gasoline or other solvents is prohibited.**
These agents may damage the skin or corrode the product.
4. **Storage and operating conditions**
Do not expose the product to strong sunlight or moisture. Otherwise, it may cause damage to the product.

SAFETY MEASURES

The device is a professional equipment for use in beauty parlors. Operation of the device only by qualified cosmetic personnel avoids undesirable consequences during treatment.

Do not disassemble the device or attempt to use it for purposes other than those indicated in this manual. All servicing and repairs must be performed only by an authorized service technician (ACTIVESHOP Ltd. company).

Do not install or operate the device near liquids, wet places or with wet hands - so that no liquids get into the device.

Immediately unplug the device from the outlet and contact the authorized service provided by ACTIVESHOP Ltd. in the following cases:

(a) liquid has entered the device, (b) the device emits a strange odor, smoke or makes strange sounds,
c) the power cord has been damaged, d) the device has been overturned or its housing has been damaged.

Do not lay any objects on the cables. Do not install the device where someone may step on the power cord.

For safety purposes, unplug the power cord from the outlet each time you finish working with the device.

Do not insert anything into the opening of the device. If a foreign object gets into the hole, remove the cord from the outlet and contact the distributor or manufacturer.

Do not place the device on an unstable cart, shelf or vehicle. If dropped, the device may be damaged.

In order to guarantee the established useful life of the device, it is recommended to keep an interval of one minute before pressing the power button again (between switching on and off and vice versa).

The power supply must comply with the specifications on the nameplate. Otherwise, there is a risk of incorrect operation of the device its damage and even burning.

CAUTION:

In case of damage to the power cord or other malfunctions, contact an authorized service center - this will avoid danger.


This appliance can be used by children at least 8 years of age and by persons with reduced physical and mental capabilities and persons with inexperience and lack of familiarity with the equipment, if supervision or instruction on the use of the equipment in a safe manner is provided, so that the associated risks are understood.

The device is not intended for play by children. Do not allow children to play with the device.

Do not leave the appliance unattended when children or persons without sufficient knowledge in the use of the appliance are near it.

This electrical appliance is not a toy. Improper use or contact with the electrical appliance may pose a health risk or lead to death, so use and store the appliance out of the reach of children.

DISPOSAL OF THE DEVICE

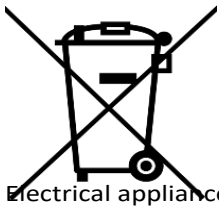

Dispose of used electrical appliances, fixtures and packaging to raw material reprocessing centers in accordance with environmental protection principles. Do not dispose of electrical appliances in household waste containers. Follow the rules for discarding this type of equipment. In order to protect the environment in accordance with Directive 2012/19/EC on waste electrical and electronic equipment, used

Electrical appliances should be collected separately, and then handed over to collection points for

collection of waste electronic and electrical equipment. At the end of their useful life, as well as old appliances contain full-value recyclable components, so they should be handed over for reprocessing. Therefore, they should be handed over to appropriate collection points, as well as they can also be handed over to ACTIVESHOP Ltd.

STORAGE AND TRANSPORT

The device should be transported at a temperature in the range of (-) 20 °C to (+) 50 °C only in its original packaging. Store the device in a well-ventilated and relative humidity of no more than 75%. Operating temperature from (+) 10 °C to (+) 40 °C. During storage and transportation processes, the device must not be thrown or dropped, exposed to mechanical shocks or atmospheric agents.

