ACTIVESHOP®

Instructions for use

Azzurro 713A hydraulic podiatry chair




Scan for different languages

CE


Importer:
ACTIVESHOP Sp. z o.o.
Ul. Graniczna 8b, bud. DC2A 54-610 Wrocław, Poland
www.activeshop.com.pl
Rev. A/06/2025

Dimensions of the chair:

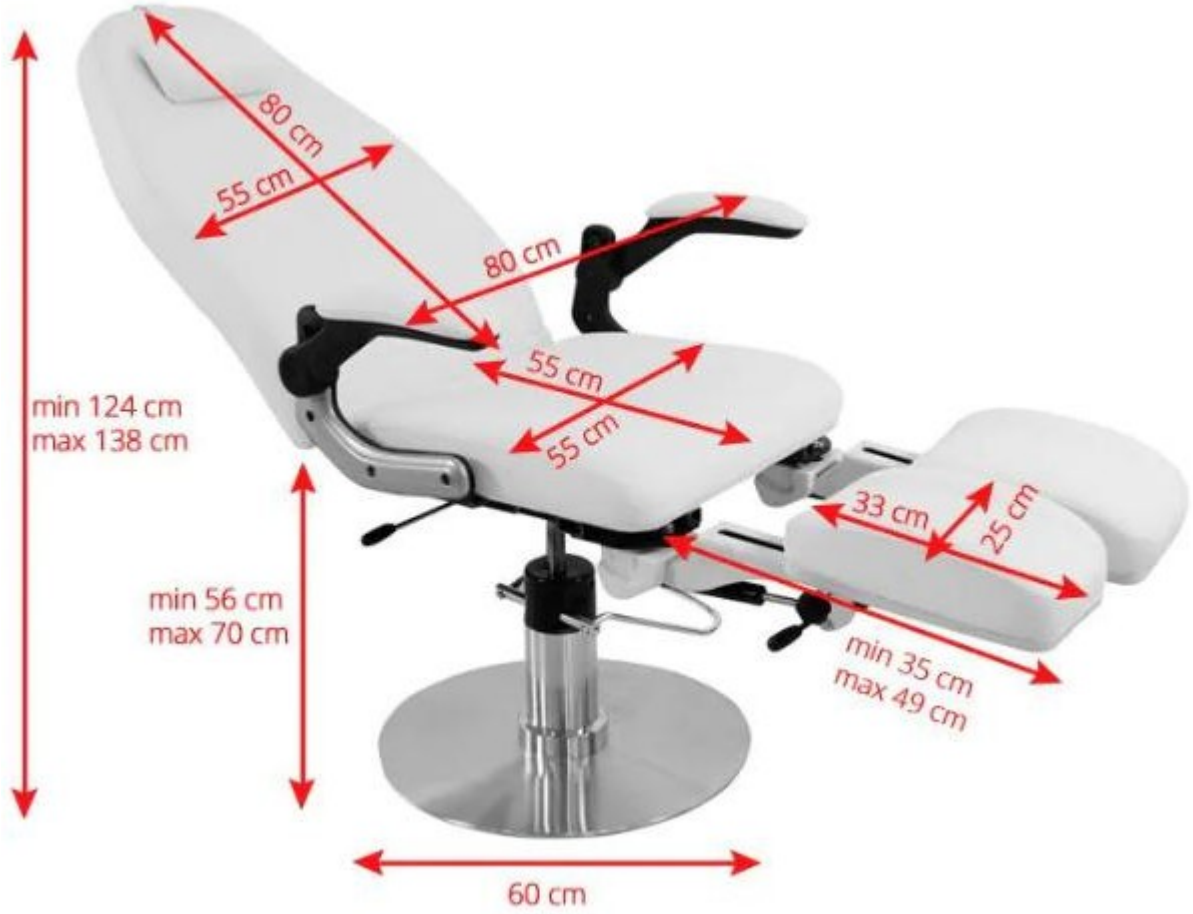


Installation:


1. Raise the actuator lever to the highest position, and then settle the seat on the actuator stem as shown.

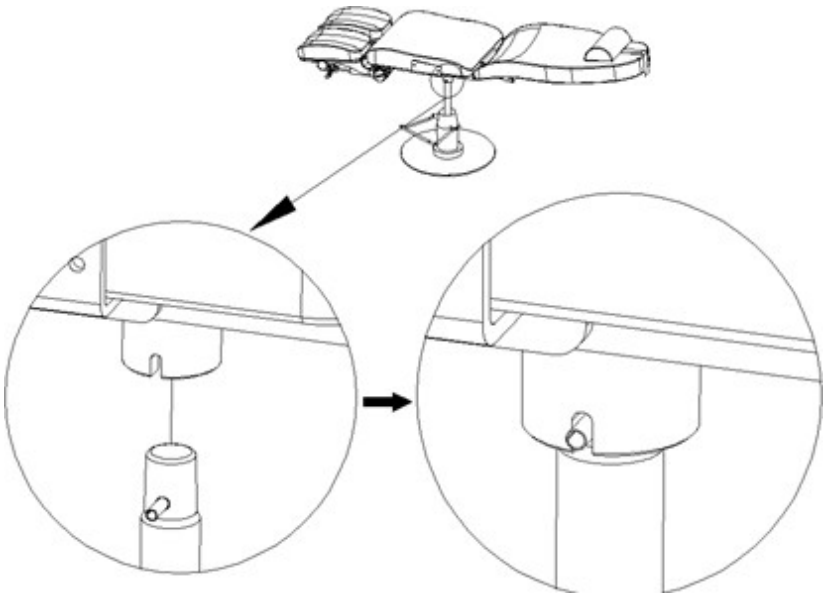

2. Align the backrest bar with the actuator hole. Insert the bolt as shown in the illustration, and then secure it with the clamp.

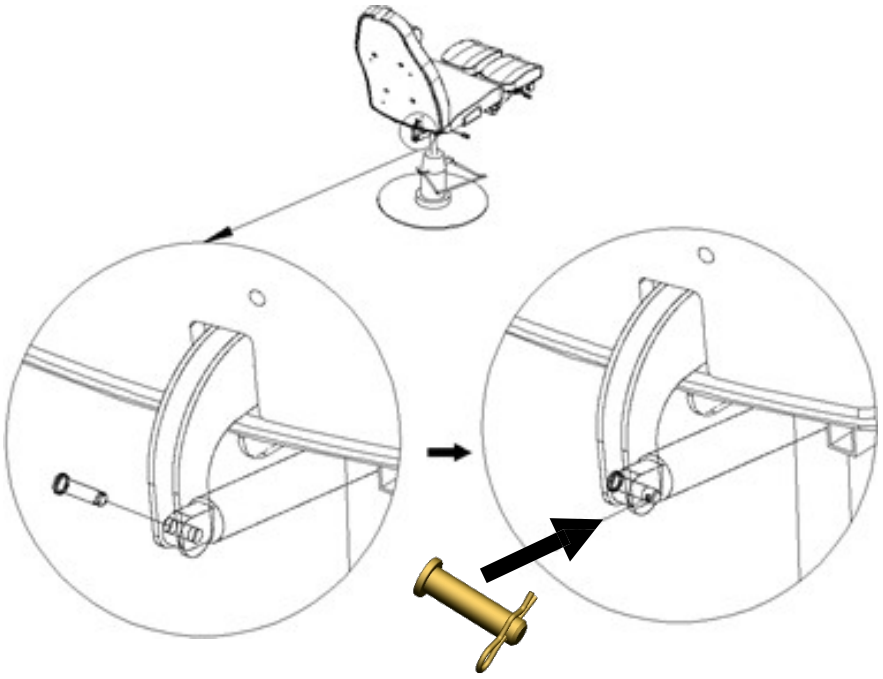


3. Mount the armrests, using the screws as shown in the illustration.

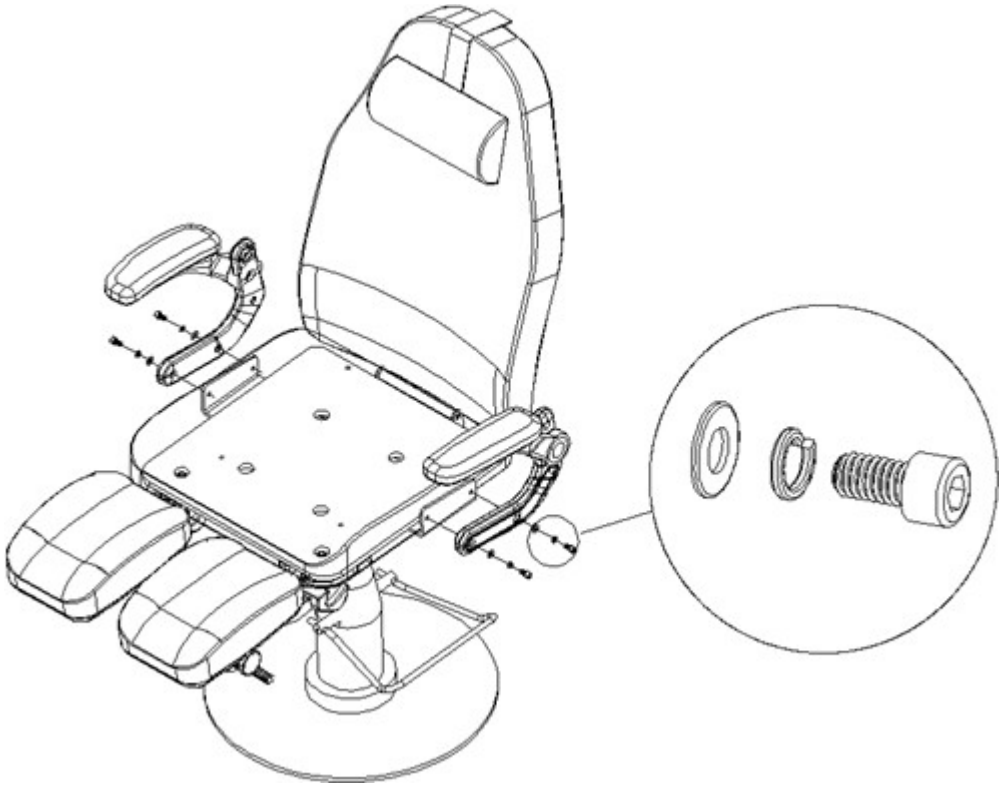


Features:


1. The seat height can be adjusted using the lever. Raising the lever causes Locking the rotation of the seat, pressing the lever lowers it.

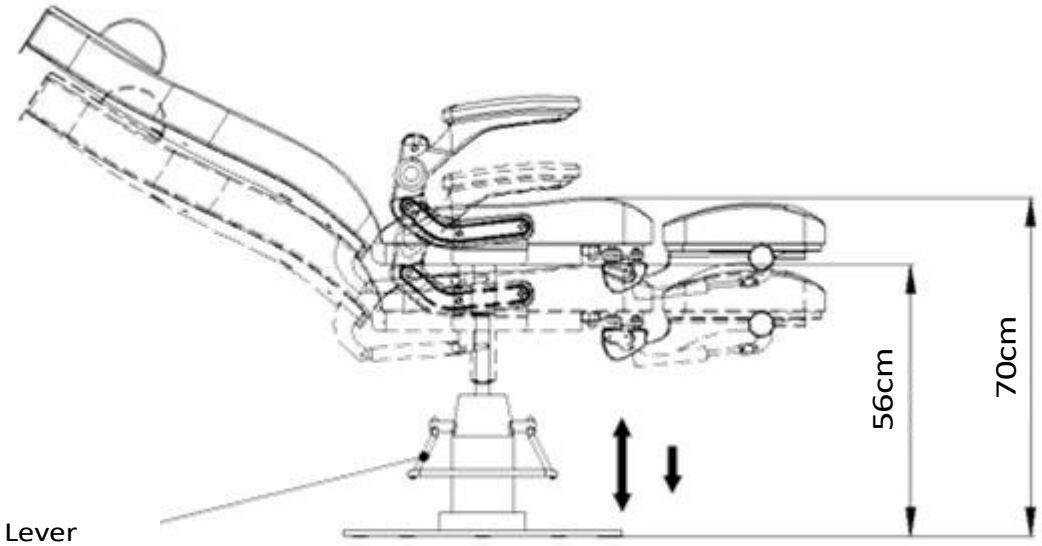

2. To tilt the backrest, press the lever and lean back, and the backrest will tilt to a convenient angle.
   The maximum tilt of the backrest is 45 degrees.

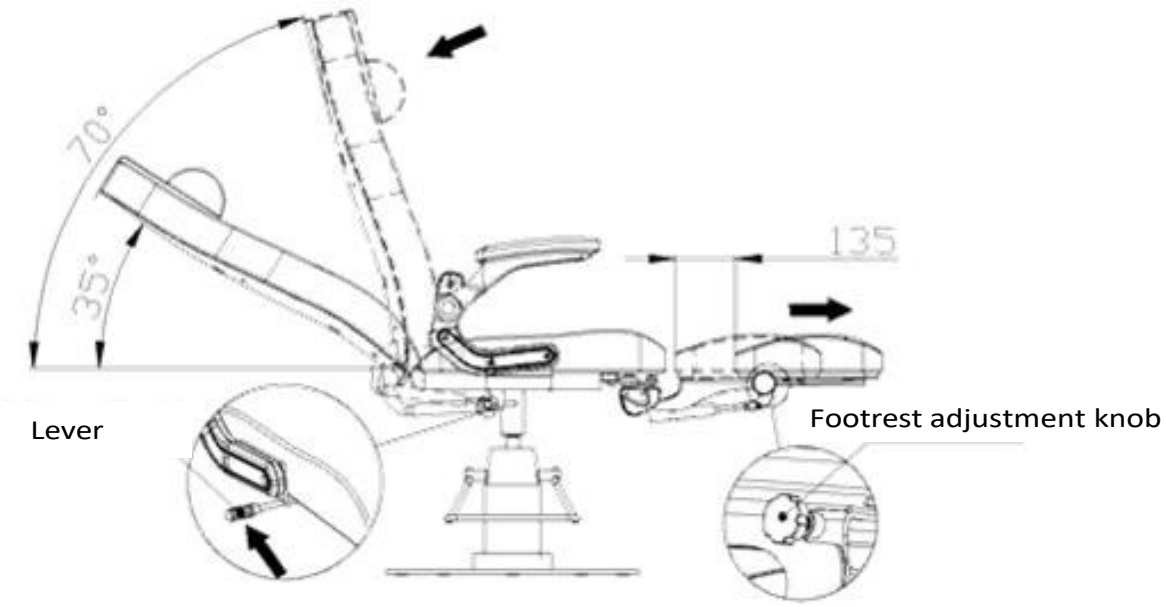


3. The footrest can be locked in 30°, 60° or 90° position.

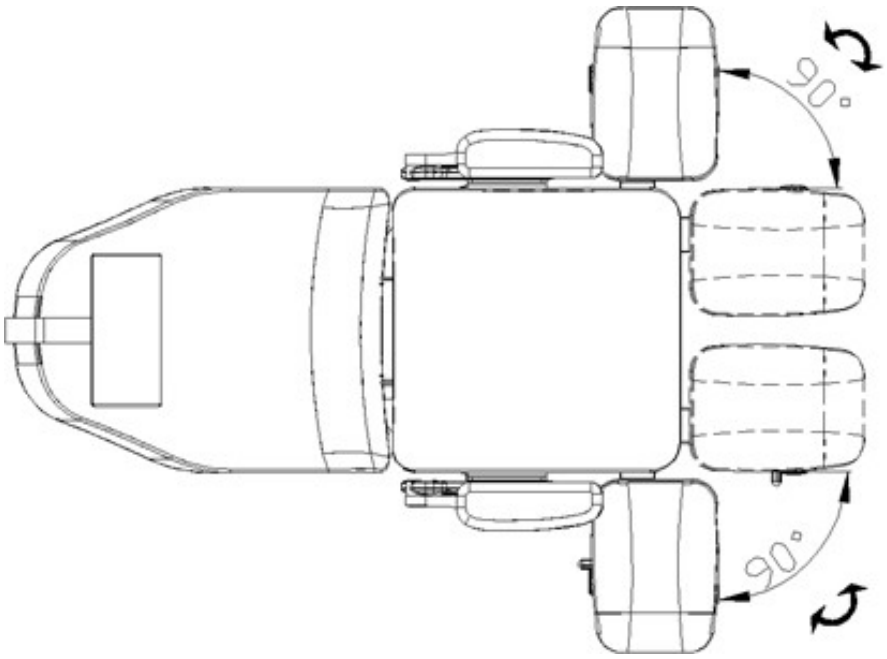

4. (a) Use the knob and simultaneously tighten the footrest to tilt it to the desired angle. Tilt range: up to 90°.

   (b) The armrest can be adjusted from horizontal to vertical position.

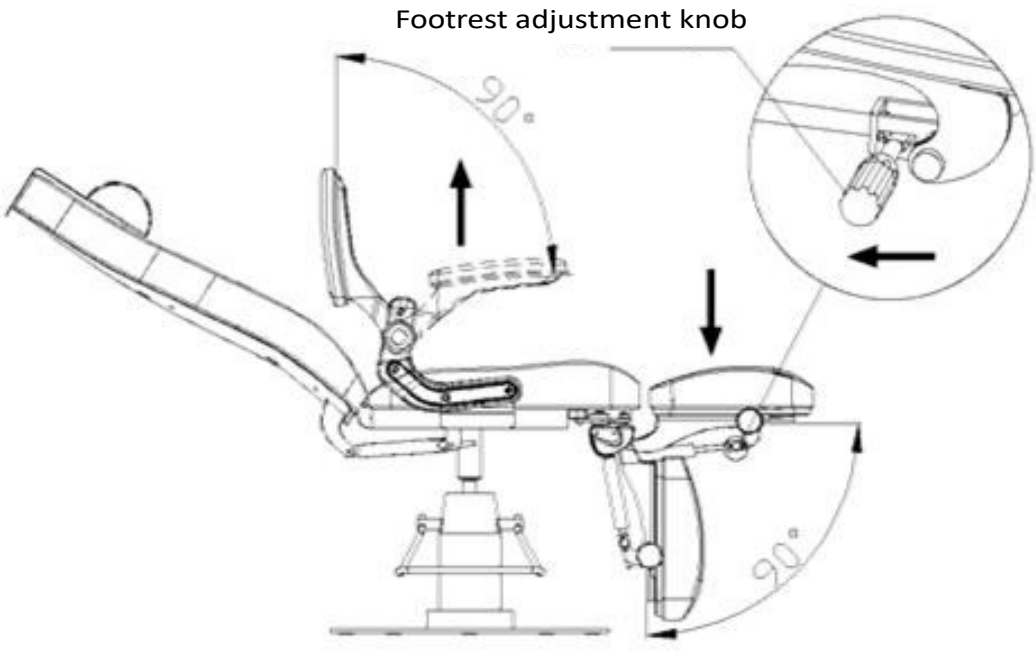

Attention.

- Ensure that children do not play near the chair without adult supervision.
- Do not put your hand or foot in the crevices of the chair.
- Do not expose the chair to direct sunlight.
- Do not push or pull the chair on an uneven floor.
- Do not stand on the chair.
- The chair should not be used by people with physical limitations, sensory impairments or mental problems.

Cleaning and maintenance:

- To keep the chair's surface in good condition, clean it regularly. Remember to clean and disinfect the chair after each use. It is recommended to perform a thorough disinfection once a month.

I Metal and plastic surfaces

- To clean metal and plastic surfaces, wipe them with a damp cloth. To clean hard-to-reach areas or corners, a small brush will work best.
- Do not use too much liquid to clean the device. Then disinfect the device with a suitable liquid.

II Stainless steel surfaces

- It will work best is cleaning with a mild detergent solution. Ammonia and most solvents can be used to remove tough stains. Avoid chlorine-based solutions.

III Painted or chrome-plated metal surfaces

- Do not use water jets (showers, high-pressure water cannons) for cleaning.
- Do not clean in high temperature and high humidity. Do not expose the device to conditions that may lead to liquid contamination.
- Do not use solvents and acids for cleaning.

IV Seat surface

- Wipe the surface of the seat evenly with a special leather cleaner, then wipe it with a clean, soft cloth.

Maintenance

1. Repairs should be performed by experienced and manufacturer-authorized personnel. Repairs by unauthorized personnel may cause injury or damage to the device, for which the manufacturer is not responsible.
2. Use only original spare parts approved by the manufacturer.
3. All maintenance and repair work must be documented.

If you find any defects, stop using the equipment and contact an authorized service center. Only experienced or authorized personnel can open or replace the actuator / control unit.

Technical data:

- Maximum load 150kg
- Weight 56.5kg
- PU/PVC leather

SAFETY MEASURES

1. The device is a professional equipment for use in beauty parlors. Operation of the device only by qualified cosmetic personnel avoids undesirable consequences during treatment.
2. Do not disassemble the device or attempt to use it for purposes other than those indicated in this manual. All servicing and repairs must be performed only by an authorized service technician (ACTIV Tomasz Pacholczyk company).
3. Do not install or operate the device near liquids, wet places or with wet hands - so,
so that no liquids get into the device.
4. Immediately remove the plug from the socket and contact an authorized service provided by ACTIV Tomasz Pacholczyk in the following cases: a) liquid has gotten into the device, b) the device emits a strange smell, smoke or makes strange sounds, c) the power cord has been damaged, d) the device has been knocked over or its housing has been damaged.
5. Do not lay any objects on the cables. Do not install the device where someone may step on the power cord.
6. For safety purposes, unplug the power cord from the outlet each time you finish working with the device.
7. Do not insert anything into the opening of the device. If a foreign object gets into the hole, remove the cord from the outlet and contact the distributor or manufacturer.
8. Do not place the unit on an unstable cart, shelf or vehicle. If dropped, the device may
be damaged.
9. In order to guarantee the established useful life of the device, it is recommended to keep an interval of one minute before pressing the power button again (between switching on and off and vice versa).
10. The power supply must comply with the specifications on the nameplate. Otherwise, there is a risk of incorrect operation of the device its damage and even burning.

CAUTION:

1. In case of damage to the power cord or other malfunctions, contact an authorized
service center - this will help avoid danger.
2. This equipment may be used by children at least 8 years of age and by persons of diminished physical and mental capabilities and persons with inexperience and lack of familiarity with the equipment, if supervision or instruction on how to use the equipment in a safe manner is provided so that the associated risks are understood.
3. The equipment is not intended for play by children. Do not allow children to play with the equipment.
4. Do not leave the appliance unattended when children or persons without sufficient knowledge of the use of the appliance are in the vicinity.
5. This electrical appliance is not a toy. Improper use or contact with the electrical appliance may pose a health risk or lead to death, so use and store the appliance out of the reach of children.

DISPOSAL OF THE DEVICE



Dispose of used electrical appliances, fixtures and packaging to raw material reprocessing centers in accordance with environmental protection principles. Electrical appliances should not be thrown into household waste containers. Follow the rules for discarding this type of equipment. In order to protect the environment in accordance with Directive 2012/19/EC on waste electrical and electronic equipment, waste electrical equipment should be collected separately and then transfer to collection points for used electronic and electrical equipment.
At the end of their useful life, as well as old appliances contain full-value recyclable components, so they should be handed over for reprocessing. Therefore, they should be handed over to appropriate collection points, as well as they can also be handed over to the ACTIV Tomasz Pacholczyk Company.

STORAGE AND TRANSPORT

The device should be transported at a temperature in the range of (-) 20 °C to (+) 50 °C only in its original packaging. Store the device in a well-ventilated and relative humidity of no more than 75%. Operating temperature from (+) 10 °C to (+) 40 °C. During storage and transportation processes, the device must not be thrown or dropped, exposed to mechanical shocks or atmospheric agents.