ACTIVESHOP®

Instructions for use

Azzurro Lux 3m beauty chair



Importer:
ACTIVESHOP Sp. z o.o.
Ul. Graniczna 8b, bud. DC2A
54-610 Wrocław, Poland
www.activeshop.com.pl
Rev. A/07/2025

Scan for different languages

Warning!

1. The footrest is not designed to support the full weight of an adult.
2. Do not sit directly on the armrest or footrest.
3. Do not lift more weight than the chair can support.
4. For safety reasons, always connect the power cord to a grounded outlet.
5. Do not tie the power cord to the device. Lifting or pulling on the power cord may damage it. Ensure that the power cord can be easily removed in an emergency.
6. Before moving the electric beauty chair, remove the power cord so that it does not get stuck between the parts of the chair.
7. If the power cord is damaged, unplug it immediately, do not use the chair, and contact an authorized service center. Use only the original power cord.
8. Do not allow anyone who is not performing the treatment to operate the remote control.
9. Do not place any objects under the chair.
10. Ensure there is sufficient space around the chair to allow it to be folded and unfolded freely.
11. Make sure there are no sharp or harmful objects around the chair that could damage the device.
12. Do not modify the device yourself. If the chair does not work properly, stop using it and contact an authorized service center.
13. Do not throw heavy objects on the chair or jump on it.

Assembly

1. Insert the headrest into the hole at the top of the backrest (Fig. 1).
2. Position and connect the left and right armrests to the pin in the chair, then rotate them 90° counterclockwise (Fig. 2).


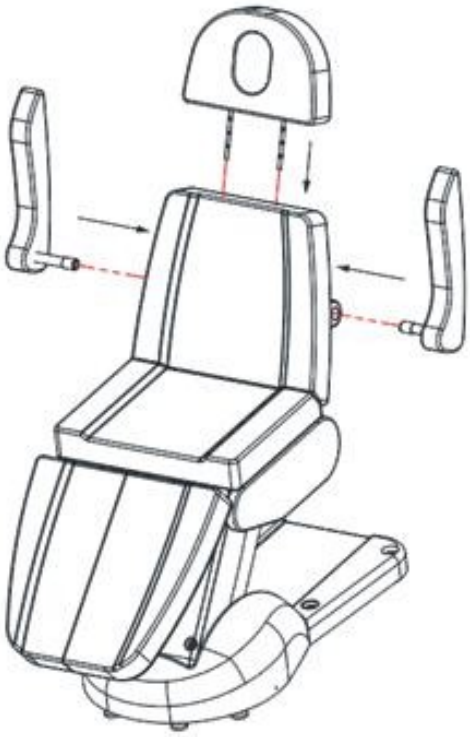


Fig. 1



Fig. 2


Control remote

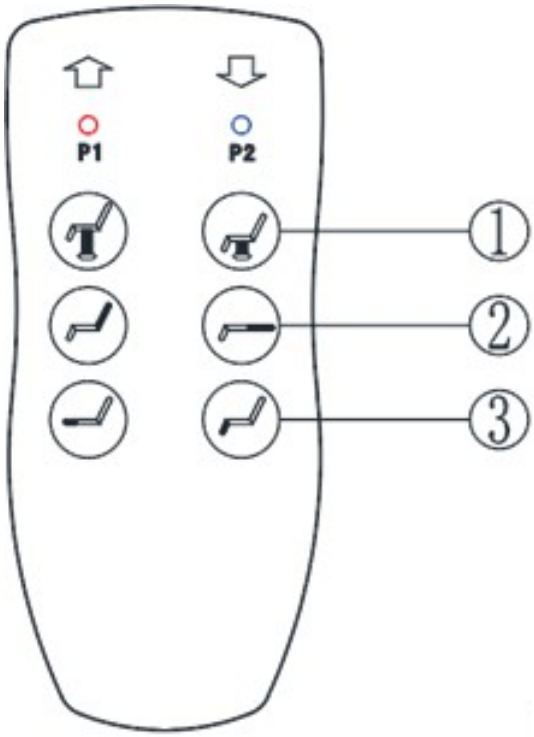


1. Press the button to raise the seat.
1. Press the button to lower the seat.
2. Press the button to increase the angle of the backrest.
2. Press the button to decrease the angle of the backrest.
3. Press the footrest tilt button to increase the angle of the footrest.
3. Press the button to decrease the angle of the footrest.

Additional features

- Pull the handle backward and hold it to adjust the angle of the footrest.

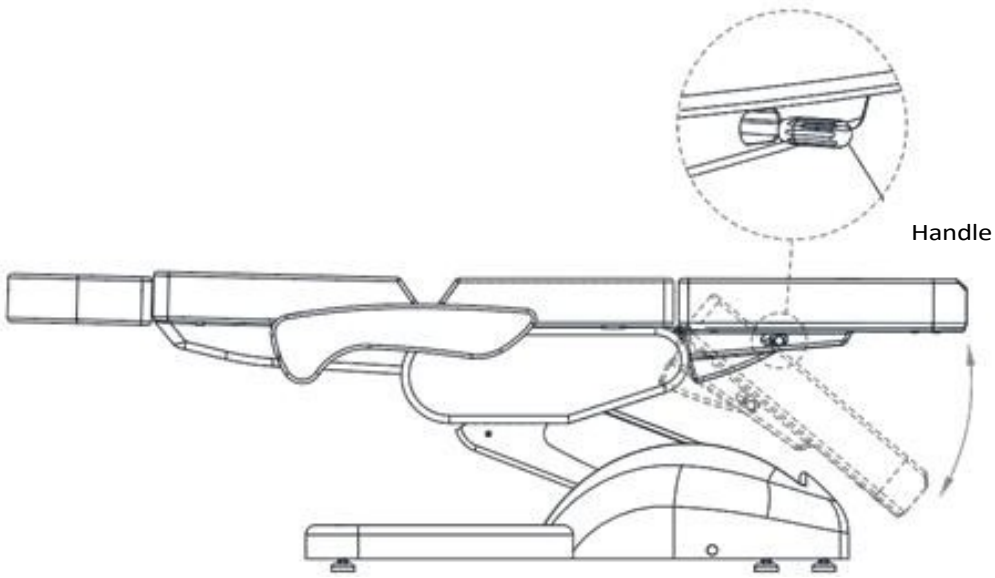

- The height of the headrest can be adjusted by pulling it up or pushing it down.

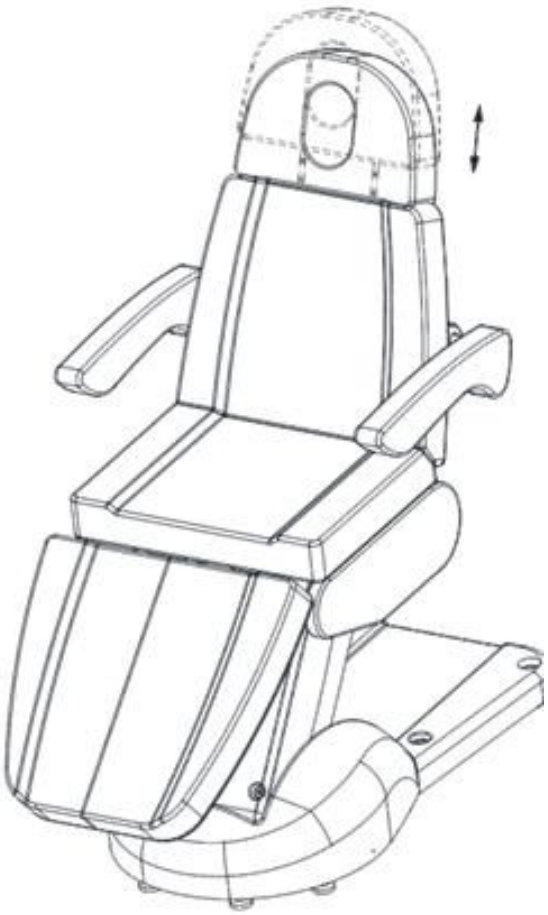

Maintenance:

Cleaning and disinfecting the device:

- Unplug the device before cleaning.
- To keep the surface of the chair in good condition, clean it regularly. Remember to clean and disinfect the chair after each use. It is recommended to perform thorough disinfection once a month.

I Metal and plastic surfaces

- To clean metal and plastic surfaces and controls, wipe them with a damp cloth.
with a cloth. A small brush is useful for cleaning hard-to-reach areas or corners. Do not use too much cleaning fluid on the device. Then disinfect the device with a suitable liquid.

II Stainless steel surfaces

- A mild detergent solution works best. Ammonia and most solvents can be used to remove stubborn stains. Avoid chlorine-based solutions.

III Painted or chrome-plated metal surfaces

- Do not use water jets (shower, high-pressure water cannons) for cleaning.
- Do not clean in high temperatures and high humidity, do not expose the device to conditions that may lead to liquid spillage.
- Do not use solvents or acids for cleaning.

IV Seat surface

- Wipe the surface of the chair evenly with a special leather cleaner, then wipe it with a clean, soft cloth.

Before servicing the device, disconnect it from the power supply.

- Repairs should be performed by experienced personnel authorized by the manufacturer.
- Repairs by unauthorized personnel may result in personal injury or damage to the device, for which the manufacturer is not responsible.
- Use only original replacement parts approved by the manufacturer.
- Do not use the electric beauty chair or accessories if they are not working properly.
- All maintenance and repair work must be documented.

Warning!

- If you find any defects, stop using the equipment and contact an authorized service center. Only experienced or authorized personnel may open or replace the actuator / control unit.
- If any part of the device is damaged, unplug the power cord and stop using it. Then notify the service center.

Troubleshooting:

Problem	Diagnosis	Solution
Buttons do not respond	- Check that the power plug is connected - Check that all transformer connection plugs are properly connected	- Unplug the power plug and plug it back in - Disconnect all connecting plugs and reconnect them
Noise and vibration occur when lifting the chair	- Loose screws	- Tighten the screws
The indicator light is not on	- The remote control cable is not properly connected or the remote control is damaged	- Reconnect the cable; if it still does not work, replace the remote control

Technical data:

Dimensions: 2000 x 810 x 610–910
mm Maximum load: 140 kg Net
weight: 76 kg
Adjustable height range: 610–910 mm
Input voltage: AC 100–130 V / 220–240 V, 50/60 Hz Output
voltage: DC 24 V

SAFETY MEASURES

1. The device is professional equipment for use in beauty salons. Only qualified beauty personnel should operate the device to avoid undesirable effects during treatment.

2. Do not disassemble the device or attempt to use it for purposes other than those specified in this user manual. All service and repair work must be performed exclusively by an authorized service technician (ACTIV Tomasz Pacholczyk).

3. Do not install or operate the device near liquids, wet areas, or with wet hands, so that no liquids can get into the device.

4. Immediately unplug the device and contact the authorized service provider ACTIV Tomasz Pacholczyk in the following cases: a) liquid has entered the device, b) the device emits a strange odor, smoke, or strange noises, c) the power cord has been damaged, d) the device has been knocked over or its casing has been damaged.

5. Do not place any objects on the cables. Do not install the device in places where someone could step on the power cord.

6. For safety reasons, unplug the power cord from the outlet after each use of the device.

7. Do not insert anything into the opening in the device. If a foreign object gets into the opening, unplug the cord from the outlet and contact the distributor or manufacturer.

8. Do not place the device on an unstable cart, shelf, or vehicle. The device may be damaged if it falls.

9. In order to guarantee the specified service life of the device, it is recommended to wait one minute before pressing the power button again (between switching on and off and vice versa).

10. The power supply must comply with the specifications on the nameplate. Otherwise, there is a risk of malfunction, damage, or even burning of the device.

CAUTION:

1. In the event of damage to the power cord or other malfunctions, contact an authorized service center to avoid danger.

2. This appliance can be used by children aged 8 years and above and by persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way so that the hazards are understood.

3. The device is not intended for use by children. Children must not be allowed to play with the device.

4. Do not leave the device unattended when children or persons who do not have sufficient knowledge of how to use the device are in its vicinity.

5. This electrical device is not a toy. Incorrect use or contact with the electrical device may pose a risk to health or lead to death, so the device must be used and stored out of the reach of children.



DISPOSAL OF THE DEVICE

Used electrical appliances, equipment, and packaging should be taken to recycling centers in accordance with environmental protection regulations. Electrical appliances should not be disposed of in household waste containers. Follow the rules for disposing of this type of equipment. In order to protect the environment natural in accordance with Directive 2012/19/EC on waste electrical and electronic equipment

Waste electrical and electronic equipment must be collected separately and then transferred to collection points for waste electrical and electronic equipment.

At the end of their useful life, old devices contain valuable recyclable components, which is why they should be sent for recycling. They should therefore be taken to appropriate collection points, or they can also be sent to ACTIV Tomasz Pacholczyk.

STORAGE AND TRANSPORT

The device should be transported at temperatures ranging from (–) 20 °C to (+) 50 °C exclusively in its original packaging. The device should be stored in a well-ventilated area with relative humidity not exceeding 75%. Operating temperature range: (+) 10 °C to (+) 40 °C. During storage and transport, do not throw or drop the device, expose it to mechanical shocks or atmospheric factors.