MOPEDIA®

ITALIANO

SOLLEVATORI KOMPASS

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MA RI800_801_804_805_807_810_814_815_820 04 B_ITA_03 2021

INDICE

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CE I Class Medical Device
REGULATION (EU) 2017/745 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL
of 5 April 2017 concerning medical devices

1. CODES

HYDRAULIC LIFTS
RI800 - HYDRAULIC actuator with payload 135 Kg
RI801 - HYDRAULIC actuator with payload 150 Kg

ELECTRICAL LIFTS
RI804 - TiMOTION actuator with payload 135 Kg
RI805 - TiMOTION actuator with payload 150 Kg
RI807 - LINAK Jumbo Home actuator with payload 135 Kg
RI810 - TiMOTION actuator with payload 180 Kg
RI814 - LINAK actuator with payload 135 Kg
RI815 - LINAK actuator with payload 150 Kg
RI820 - LINAK actuator with payload 180 Kg

2. INTRODUCTION

Thank you for choosing a patient lift/stand up lift from the **KOMPASS by MOPEDIA** product range by MORETTI S.p.A.. Their design and quality are a guarantee of comfort, safety and reliability. **KOMPASS by MOPEDIA** patient lifts/stand up lifts by MORETTI S.p.A. are designed and built to meet all your demands for a practical, correct use. This user manual provides for some suggestions as to how correctly use the lift you have chosen and offers a lots of valuable advice for your safety. Please read carefully through the manual before using the patient lift/stand up lift. Should you have any queries, please contact your retailer for advice and assistance.

3. INTENDED USE

KOMPASS by MOPEDIA patient lifts/stand up lifts are intended to lift and move people with motor difficulties. An assistant operates the machine. The lifts cannot be used by only one person. MORETTI S.p.A. refuses all liability for damages caused by improper use of the device or by its use for a purpose other than that one specified in this manual.

4. DECLARATION OF CONFORMITY

Moretti S.p.A. declares under its sole responsibility that the product made and traded by Moretti S.p.A.and belonging to the group of PATIENT LIFTS MOPEDIA complies with the provisions of the regulation 2017/745 on MEDICAL DEVICES of 5 April 2017.
For this purpose, Moretti S.p.A. guarantees and declares under its sole responsibility what follows:

1. The devices satisfy the requirements of general safety and performance requested by the Annex I of regulation 2017/745 as laid down by the Annex IV of the above mentioned regulation.
2. The devices ARE NOT MEASURING INSTRUMENTS.
3. The devices ARE NOT MADE FOR CLINICAL TESTS.
4. The devices are packed in NON-STERILE BOX.
5. The devices belong to class I in accordance with the provisions of Annex VIII of the above mentioned regulation
6. Moretti S.p.A. provides to the Competent Authorities the technical documentation to prove the conformity to the 2017/745 regulation, for at least 10 years from the last lot production.

ENGLISH

Note: Complete product codes, the manufacturer registration code (SRN), the UDI-DI code and any references to used regulations are included in the EU declaration of conformity that Moretti S.p.A. releases and makes available through its channels.

5. GENERAL WARNINGS

WARNING!
Always pay close attention to the presence of moving parts that could cause limb entrapment and injury.

- For proper use of the device refer to this manual, checking with a particular attention at the warnings before use.
- For proper use of the device to consult your doctor or therapist.
- Keep the packaged product away from any source of heat because the packaging is made of cardboard.
- The useful life of the accessory is determined by wear of parts that cannot be repaired and/or replaced.
- Be careful when children are nearby.
- Do not exceed the maximum payload.
- The user and/or the patient will have to report any serious accident that have occurred related the device to the manufacturer and appropriate authority of the State which the user and/or patient belongs to.

6. SYMBOLS

Product code

Unique Device Identification

CE mark

Manufacturer

Batch Lot

Read the instruction manual

Medical Device

Conditions of disposal

Attention, should read the instructions

Production dates


Waste Electrical & Electronic Equipment (WEEE)


IP54 Protection against dust and water splashes

B Type


7. DIMENSIONS


7.1 Dimensions of structure with payload 135 Kg
RI800 - RI804 - RI807 - RI814

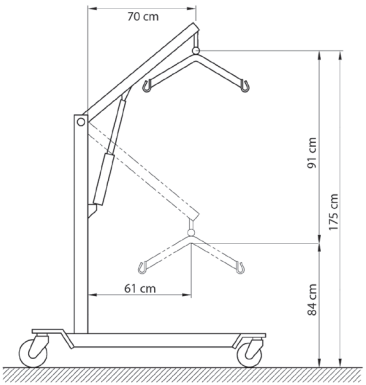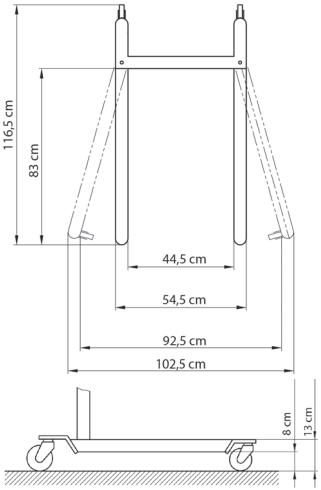
7.2 Dimensions of structure with payload 150/180 Kg
RI801 - RI805 - RI810 - RI815 - RI820

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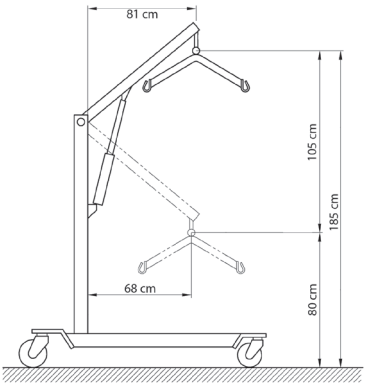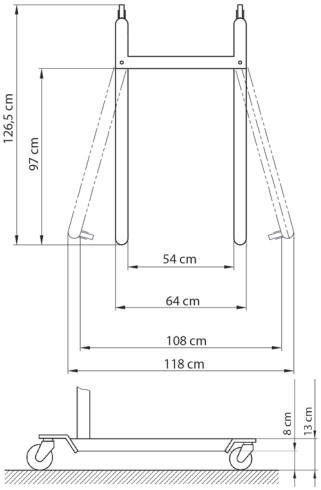
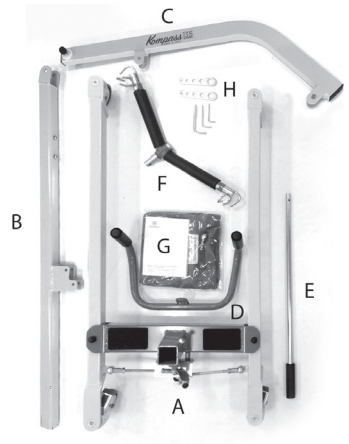
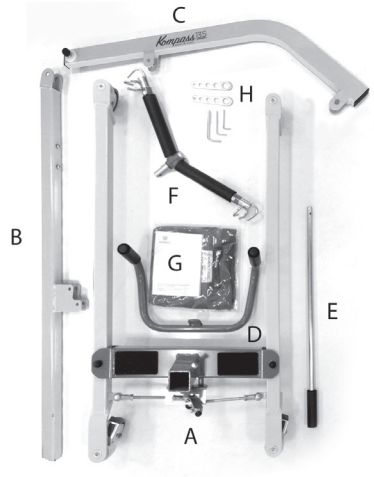
9. DESCRIPTION OF THE ACTIVATION KITS

9.1 Hydraulic actuator (135/150 KG)
Suitable for the models RI800 - RI801

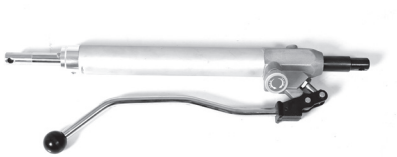
9.2 TiMOTION electrical kit (135/150 Kg)
Suitable for the models **RI804 - RI805**

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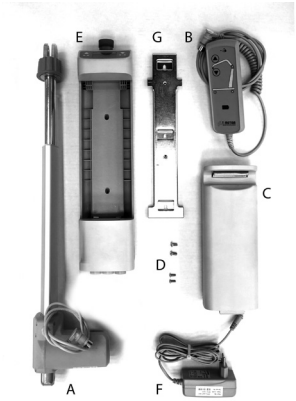
9.3 TiMOTION electrical kit (180 Kg)
Suitable for the model **RI810**

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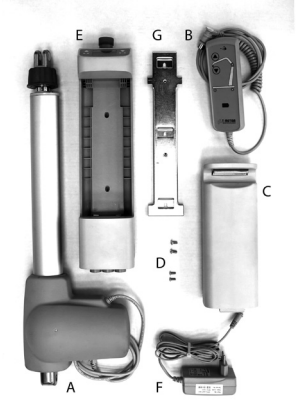
9.4 LINAK Jumbo Home electrical kit (135 Kg)
Suitable for the model **RI807**

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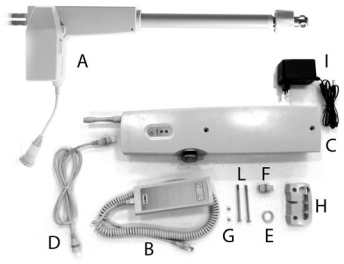ENGLISH


9.5 LINAK electrical kit (135/150/180 Kg)
Suitable for the models **RI814 - RI815 - RI820**

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
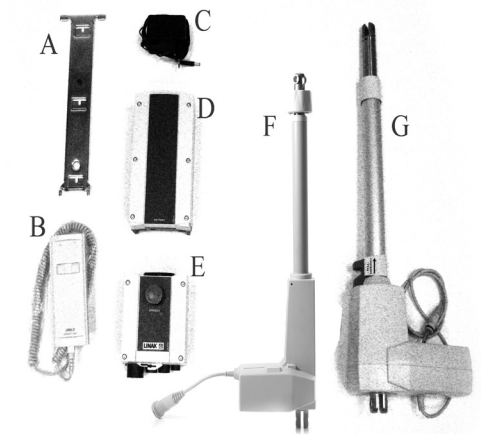
10. ASSEMBLY INSTRUCTIONS FOR LIFT STRUCTURES (FOR ALL THE MODELS)


1. Unscrew the screw at the base of the mast.
2. Fit the mast into the provided cavity of the base. (Fig.4)
3. Fit the screw that blocks the mast into the base and screw it well. (Fig.4)
4. Connect the boom to the mast with the provided screws. (Fig.5a-5b)
5. Fit the handle into the mast and block it with the 2 screws. (Fig.6)
6. Fit the lever for opening the base and block it with its screw. (Fig.7)
7. Fit the lifting beam into the corresponding hook, placed at the end of the boom. Fit the plastic washers on both sides of the anchor bolt. (Fig.8)
8. The structure will be completely assembled. Install the hydraulic/electrical actuator suitable for the chosen patient lift

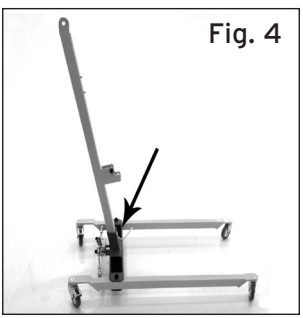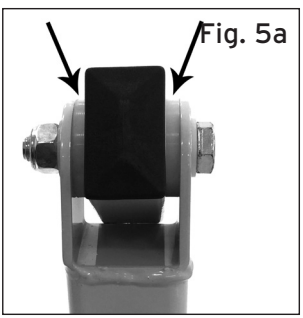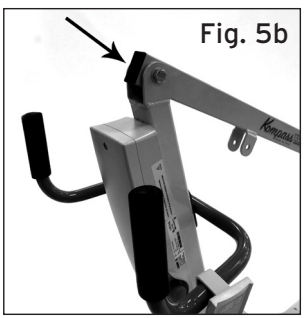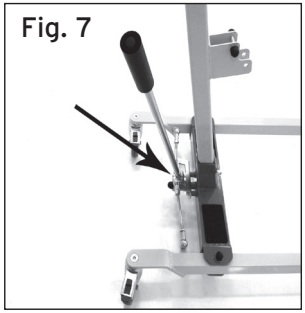
11. ASSEMBLY INSTRUCTIONS ACTUATORS


11.1 Assembly HYDRAULIC actuator (suitable for RI800 - RI801)



1. Place the hydraulic actuator as in the picture, blocking well its upper and lower part by auto-blocking screws and bolts. (Fig.9)


11.2 Assembly TiMOTION electrical kit (suitable for RI804 - RI805 - RI810)

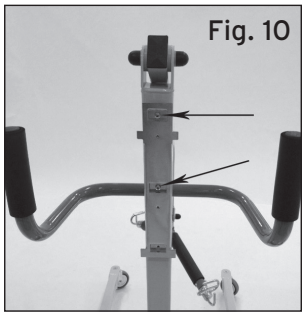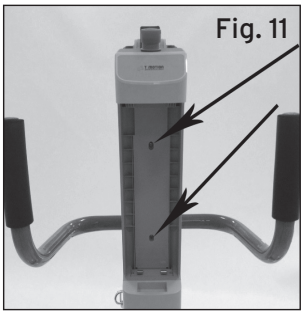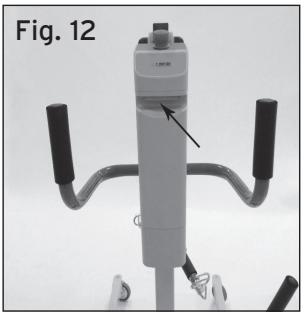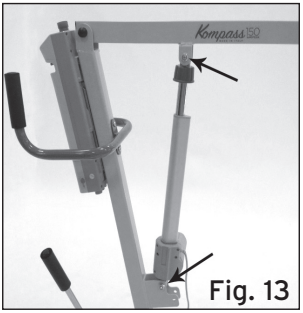
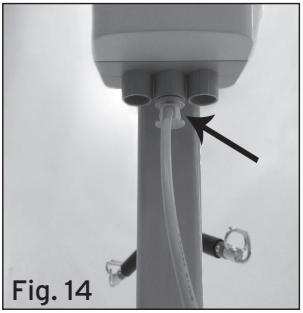
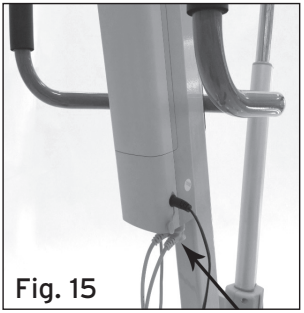
1. Secure the support mount of the control unit to the mast with the provided screws, using the fitted cavities. (Fig.10)
2. As soon as you have assembled the mount, fit the control unit on it, letting it slide from top to bottom and blocking it with the screws. (Fig.11)
3. Fit the battery into the control unit. (Fig.12)
4. Remove the two actuator mounting bolts from the mast boom. Fit the actuator into it and secure it with the two bolts previously removed. Cover the bolt heads with the provided plastic cap. (Fig.13)
5. Wire up the actuator connector to the central socket of the control unit and the remote control to the right socket. (Fig.14-15)


11.3 Assembly LINAK Jumbo Home electrical kit (suitable for RI807)

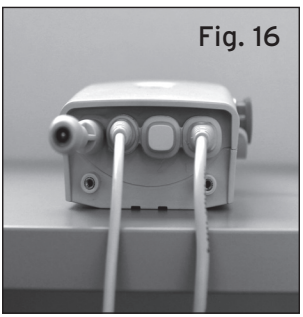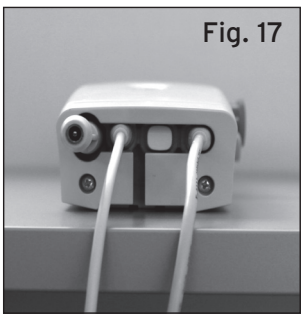1. Connect the end of the extension cable into the right part of the unit. Wire up the plug of the remote control into the correspondent socket of the control unit, making sure the guide lug seats properly. (Fig.16)
2. Fit the cap into the central connector of the control unit and then place the protection cover, using the two little screws M4x10mm. (Fig.17)
3. Remove the two screws on the mast and on the mast boom.
4. Fit the actuator and secure it well with the two screws previously removed. (Fig.18)

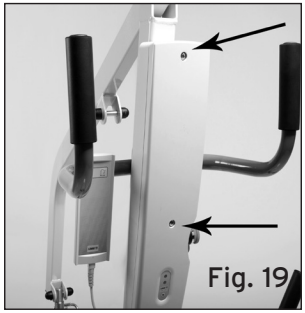Fig. 19

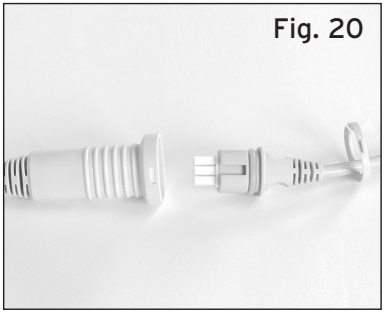
Fig. 20

Fig. 21

5. Fit the complete kit, including the battery, directly on the patient lift mast, using the two long screws M5x80mm. (Fig.19)

Note:: FIRSTLY, LIGHTLY SCREW BOTH SCREWS AND THEN TIGHTEN THEM DOWN EQUALLY UNTIL BOTH ARE FULLY SECURED

6. Connect the extension cable end to the electrical actuator, making sure the guide lug seats properly. After that, secure the connection with the plastic washer. (Fig.20-21)

11.4 Assembly LINAK electrical kit (suitable for RI814 - RI815 - RI820)

Fig. 22

Fig. 23

Fig. 24

1. Secure the control unit support mount to the mast with the provided screws, using the fitted cavities. (Fig.22)
2. Fit the battery into the control unit. (Fig.23)
3. Fit the elevation actuator between the boom and the mast and secure it firmly with screws and bolts. (Fig.24)

Fig. 25

4.Wire up the electrical connector of the actuator to the control unit socket "1" - make sure it is fully inserted. Wire up the remote control electrical connector to the large socket on the control unit - make sure the guide lug seats properly and then push it all the way in. (Fig.25)

12. INSTRUCTIONS FOR USE

12.1 Instructions for using HYDRAULIC models RI800 - RI801

1. Before lifting a patient, open the base legs to the maximum width possible, gently polling the opening lever towards you and secure it in the open position
2. Block the rear wheel brakes
3. Keep the patient centred between the base legs and always turned towards the assistant
4. Hook the cradle straps securely to the lifting beam
5. Make sure the hydraulic pump valve is fully closed (turn it in the direction marked as CLOSE)
6. Activate the pumping lever to start lifting
7. ALWAYS MONITOR THE PATIENT'S POSITION THROUGHOUT THE LIFTING PROCESS!
8. Lift the patient to the desired height, keeping him facing his assistant at all times
9. Release the rear wheels and move the patient
10. Before lowering the patient, block the rear wheels again
11. Start lowering the patient by slowly turning the hydraulic pump valve in the direction marked as OPEN

Note: To lower the patient slowly, open the valve just a little bit; opening it further always in the direction marked as "OPEN", if you want to lower the patient more rapidly

12. ALWAYS MONITOR THE PATIENT'S POSITION THROUGHOUT THE LIFTING PROCESS!

ALWAYS HELP THE PATIENT DURING THE LOWERING PROCESS SO AS TO MAINTAIN THE CORRECT SEATED POSITION!

WARNING!
While lowering the patient, please remove any barrier that can impede the downward movement. If the case, lift up the patient again and remove what causes the impediment for lowering.

12.2 Instructions for using TiMOTION electrical models RI804 - RI805 - RI810

BEFORE USING THE ELECTRICAL LIFT YOU ARE RECOMMENDED TO CHARGE THE BATTERY FOR ABOUT 6 HOURS

1. Before lifting a patient, open the base legs to the maximum width possible, gently pulling the opening lever towards you and secure it in the open position (only if there is a mechanical opening of the base legs)
2. Block the rear wheel brakes
3. Keep the patient centred between the base legs and always turned towards the assistant
4. Hook the cradle straps securely to the lifting beam
5. Press the up-facing arrow key in order to lift the patient (for the model RI810 the led display on the control unit will show the boom moving)
6. ALWAYS MONITOR THE PATIENT'S POSITION THROUGHOUT THE LIFTING PROCESS!
7. Lift the patient to the desired height, keeping him facing his assistant at all times
8. Release the rear wheels and move the patient
9. Before lowering the patient, block the rear wheels again
10. In order to lower the patient, press the down-facing arrow key (for the model RI810 the led display on the control unit will show the boom moving)
11. ALWAYS MONITOR THE PATIENT'S POSITION THROUGHOUT THE LIFTING PROCESS!

ALWAYS HELP THE PATIENT DURING THE LOWERING PROCESS SO AS TO MAINTAIN THE CORRECT SEATED POSITION!

WARNING!
While lowering the patient, please remove any barrier that can impede the downward movement. If the case, lift up the patient again and remove what causes the impediment for lowering.




EMERGENCY STOP
The control unit is equipped with a red "STOP" button which instantly deactivates all controls when pressed. To restore normal operation, turn the button clockwise and pull it outwards.

EMERGENCY KEYS
On the control unit there are two keys (one with an up-facing arrow and the other one with a down-facing arrow) that can be pressed instead of those ones on the remote control to complete all the operations when the remote control fails to work properly.

COLOR DISPLAY (Only for the model RI810)
This model is equipped with an electronic control unit with color graphic display.

LEFT KEY: Pressing the left key for more than 0,5 seconds, the boom lifting function can be activated		**RIGHT KEY:** Pressing the right key for more than 0,5 seconds, the boom lowering function can be activated

If one of the two buttons is pressed for less than 0,5 seconds, one of the following functions will be activated:

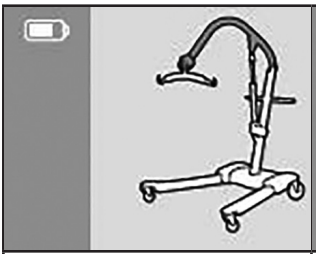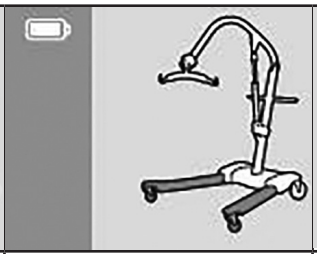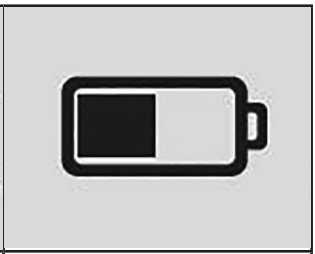		
Activation lifting boom (valid for all the models)	Electrical activation base legs opening **(only with an electrical actuator for opening base legs)**	Residual battery capacity

MAINTENANCE	MAINTENANCE
	
Months	Cycles
Time of use in MONTHS	Cycles done

NOTICE (Only for the model RI810)

ATTENTION	ATTENTION
	
Recharge Or Exchange	Overloaded
Low battery capacity Recharge or replace the battery	**Excess load alarm** The load lifted by the patient lift exceeds the maximum level allowed. Remove the load and try again.

MECHANICAL LOWERING OF THE ACTUATOR

If using the manual/mechanical lowering, please make sure that it will be done in total safety conditions for the patient. In order to activate the manual/mechanical lowering, pull downwards the red ring placed in the upper part of the actuator and turn it counter-clockwise. Once finished the lowering process, take the red ring back to its resting position.

INDICATION OF THE BATTERY STATE OF CHARGE

There is a 3-Led indicator for the battery state of charge. This indicator is placed on the remote control. 3 Leds on indicate a complete charge of the battery, whereas one just Led on indicates the battery is out of power and needs to be charged.

RECHARGING BATTERIES

- Connect the adapter to the main socket and then to the socket on the control unit, after removing the protection cap
- On the charger adapter there is a led that during the charging process becomes orange color. Charge the battery until the led becomes green color.
- Before activating the lift again, make sure to disconnect the battery charger!
- Charge the batteries for at least 6 hours before using the lift

12.3Instructions for using LINAK electrical model RI807
BEFORE USING THE ELECTRICAL LIFT YOU ARE RECOMMENDED TO CHARGE THE BATTERY FOR ABOUT 6 HOURS

1. Before lifting a patient, open the base legs to the maximum width possible, gently pulling the opening lever towards you and secure it in the open position
2. Block the rear wheel brakes
3. Keep the patient centred between the base legs and always turned towards the assistant
4. Hook the cradle straps securely to the lifting beam
5. Press the "lifting" button to lift the patient
 (if fine, a green light onto the control unit turns on)
6. ALWAYS MONITOR THE PATIENT'S POSITION THROUGHOUT THE LIFTING PROCESS!
7. Lift the patient to the desired height, keeping him facing his assistant at all times
8. Release the rear wheels and move the patient
9. Before lowering the patient, block the rear wheels again
10. In order to lower the patient, press the "lowering" button
 (if fine, a green light onto the control unit turns on)
11. ALWAYS MONITOR THE PATIENT'S POSITION THROUGHOUT THE LIFTING PROCESS!

ALWAYS HELP THE PATIENT DURING THE LOWERING PROCESS SO AS TO MAINTAIN THE CORRECT SEATED POSITION!

EMERGENCY STOP
The control unit is equipped with a red "STOP" button which instantly deactivates all controls when pressed. To restore normal operation, turn the button clockwise and pull it outwards.

EMERGENCY LOWERING
On the control unit, under the two Led lights, there is a down-facing arrow key that can be pressed for an emergency lowering when the remote control fails to work properly.

MECHANICAL LOWERING OF THE ACTUATOR
If using the manual/mechanical lowering, please make sure that it will be done in total safety conditions for the patient. In order to activate the manual/mechanical lowering, pull downwards the plastic ring placed in the upper part of the actuator and turn it counter-clockwise. Once finished the lowering process, take the plastic ring back to its resting position.

RECHARGING BATTERIES
• Connect the adapter to the main socket and then to the socket on the control unit.

Note: WHEN THE BATTERY IS CHARGING, A LED LIGHT ORANGE COLOR ONTO THE BATTERY COVER TURNS ON AND ALL THE CONTROLS ARE AUTOMATICALLY DISACTIVATED

• Before activating the lift again, make sure to disconnect the battery charger!
• Charge the batteries for at least 6 hours before using the lift

12.4 Instructions for using LINAK electrical models RI814 - RI815 - RI820
BEFORE USING THE ELECTRICAL LIFT YOU ARE RECOMMENDED TO CHARGE THE BATTERY FOR ABOUT 6 HOURS

1. Before lifting a patient, open the base legs to the maximum width possible, gently pulling the opening lever towards you and secure it in the open position
2. Block the rear wheel brakes
3. Keep the patient centred between the base legs and always turned towards the assistant

4. Hook the cradle straps securely to the lifting beam
5. Press the up-facing arrow key to lift the patient
6. ALWAYS MONITOR THE PATIENT'S POSITION THROUGHOUT THE LIFTING PROCESS!
7. Lift the patient to the desired height, keeping him facing his assistant at all times
8. Release the rear wheels and move the patient
9. Before lowering the patient, block the rear wheels again
10. In order to lower the patient, press the down-facing arrow key
11. ALWAYS MONITOR THE PATIENT'S POSITION THROUGHOUT THE LIFTING PROCESS!

ALWAYS HELP THE PATIENT DURING THE LOWERING PROCESS SO AS TO MAINTAIN THE CORRECT SEATED POSITION!

WARNING!
While lowering the patient, please remove any barrier that can impede the downward movement. If the case, lift up the patient again and remove what causes the impediment for lowering.

AUTOMATIC ANTI-CRUSHING SYSTEM
The electrical kit is equipped with an automatic anti-crushing system, which stops the lowering of the lifting boom if there is an accidental contact between the patient and any barrier.

EMERGENCY STOP
The control unit is equipped with a red "STOP" button which instantly deactivates all controls when pressed. To restore normal operation, turn the button clockwise and pull it outwards.

EMERGENCY KEYS
Underneath the red STOP button of the control box, there is a small button that can be used to low patients when the remote control fails to operate properly. This button can only be operated with the aid of a narrow tipped object like the tip of a pen.

ACTUATOR QUICK RELEASE (only model RI820)
If using the emergency lowering system, please make sure that it will be done in total safety conditions for the patient. The actuator of the lift model RI820, payload 180 Kg, has got a quick release system to activate the patient's emergency lowering when the electrical parts do not work properly, or when the patient requires immediate assistance. In order to activate the emergency lowering, please lift up the red protection cover.

Note: THE EMERGENCY QUICK RELEASE SYSTEM HAS GOT A FULLY MECHANICAL OPERATION AND IT IS CALIBRATED TO WORK WITH A MINIMUM LOAD APPLIED. SO, THE LIFTING BEAM WILL NOT BE WORKING IF IT WILL BE USED WITHOUT THE PATIENT SUPPORTED BY THE CRADLE!

RECHARGING BATTERIES
- Connect the adapter to the main socket and then to the socket on the control unit, after removing the protection cap
- On the charger adapter there is a Led light that during the charging process becomes red color. Charge the battery until the led becomes green color.
- Before activating the lift again, make sure to disconnect the battery charger!
- Charge the batteries for at least 6 hours before using the lift
- It is possible to charge the battery even using the wall charger **(OPTIONAL)**

13. WARNINGS

- The lift is a device designed for moving patients with motor problems for short distances, for instance from a bed to a wheelchair or comfortable chair, and back. It must not be used for transporting the patient.
- While lifting the patient, make sure the base legs are as far apart as possible, the control lever is securely locked and the rear wheel brakes are secured.
- While moving the patient, please avoid crossing barriers that might cause the lift to tip over.

WARNING!
Always check that the cradle is in good condition. Do not use cradles showing signs of damage or wear.

14. CLEANING

To clean the cradle and the structure use a solution of alcohol and neutral soap (in equal parts) mixed with water.

15. MAINTEINANCE

- It is extremely important to check the lift before each use.
- More specifically, please check that all screws and bolts are well tightened.
- Please, check the lift can move freely with the wheels unlocked by the brakes.
- Check the rear wheel brakes function correctly.
- It is very important to regularly check the connection pins and stop rings are correctly located (the pins connect the actuator to the upper and lower part of the lift, and the lifting beam where you hook the cradle).
- The mechanical parts must be lubricated with a drop of oil in order to facilitate their sliding and resistance to wear. In particular, pay attention to the mechanical retractor for the legs of the patient lift. Lubricate as indicated above every two or three months.

16. SPARE PARTS

For spare parts and accessories, please always refer to the MORETTI master catalogue.

ONLY USE ORIGINAL MORETTI ACCESSORIES. MORETTI S.P.A. COMPANY SHALL NOT BE LIABLE FOR ANY INCIDENT CAUSED BY IMPROPER USE.

17. CONDITIONS OF DISPOSAL 

When disposing of the device never use the normal return systems of municipal solid waste. It is recommended instead to disposal through the joint municipal ecological islands for the expected operations of recycling the materials used.

17.1 Warnings for the proper disposal of the product, as required by the European Directive 2012/19/UE 

At the end of its service life, the product must not be disposed of together with normal municipal waste. It must be delivered to the municipal sorted waste collection centre or to the resellers, in case they provide for this service. Sorted waste disposal reduces the environmental and health consequences of an improper disposal. Furthermore, the materials which make up the product may be recycled, thus resulting in energy and resource savings. In order to highlight the compulsory sorted disposal of the electro-medical equipment, the product itself bears the barred bin symbol.

17.2 Treatment of dead batteries (Directive 2006/66/CE) 

The symbol on the product indicates that the batteries must not be treated like common household waste. Make sure the batteries are properly disposed; this prevents the possible environmental and health consequences of an improper disposal. Recycling helps to preserve natural resources. Deliver dead batteries to a proper recycling centre. For more detailed information regarding the disposal of dead batteries or of the product itself, please contact the municipality, the local waste disposal centre or the reseller from whom you purchased the device.

18. WARRANTY

Every Moretti product is guaranteed for 2 (two) YEARS from the purchase date, except the restrictions specified below. The warranty shall not apply in case of damage caused by improper use, abuse, alterations of the product or any other use not in compliance with the user manual. The proper use of the product is indicated in the user manual. Moretti S.p.A. shall not be liable for any damages or personal injuries due to incorrect assembly/improper use of the device, not compliant with the user manual. Moretti does not guarantee its products against damages or faults in the following cases: natural disasters, non authorized maintenance or reparations, damages caused by electrical supply problems (if provided), use of parts or components not provided by Moretti, failure in accomplishing the guidelines and instructions for use, non authorized alterations, damages due to the transport (different from Moretti original delivery) or failure in maintaining the device as indicated in the manual. The components subject to wear are not covered by this warranty if the damage is due to the normal use of the product.

18.1 Warranty for rechargeable batteries (if provided)

The original batteries and the spare batteries are covered by a (ninety) 90-days warranty as far as their performance is concerned and by a (six) 6-months warranty as far as manufacturing faults are concerned or as required by law. Warranty on fully charged batteries left in an unused state for more than (three) 3 consecutive months is automatically void. Warranty on fully discharged batteries left in an unused state for more than (three) 3 consecutive days is automatically void.

18.2 Reparations

If a Moretti product shows material or manufacturing faults during the warranty period, Moretti S.p.A., together with the client, will consider whether the product fault is covered by warranty. Moretti S.p.A., at its unquestionable discretion, may replace or repair the item under warranty in a Moretti reseller point or at its own qualified offices. Labour costs for repairing the product may be at Moretti's expense if it is determined the reparation is under warranty. A reparation or replacement does not renew nor extend the warranty.

A product, not covered by warranty, may be returned for reparation only after receiving prior authorization from Moretti client service. Labour and delivery costs for an out-of-warranty reparation will be fully at the client's or reseller's expense. Reparations of out-of-warranty products are guaranteed for 6 (six) months, from the day the client receives the repaired product.

The client is informed if, after examining and testing a returning product, Moretti decides the product is not defective. The product will be sent back to the client and all the return delivery costs will be at his expense.

18.3 Spare parts
The original Moretti spare parts are guaranteed for 6 (six) months from the day you receive them.

18.4 Exoneration clauses
Except what clearly specified in this warranty and in accordance with the law, Moretti S.p.A. does not offer its clients any other declaration, warranty or condition, expressed or implicit, including any declaration, warranty or condition of merchantability, of non violation or non interference, of suitability for a specific aim. Moretti S.p.A. does not guarantee that the use of Moretti products will be continuous and faultless. The duration of possible implicit warranties that may be applied in accordance with the law is limited to the warranty period. Some states or countries do not allow limitations to the duration of an implicit warranty, or the exoneration, or the limitations of accidental or indirect damages related to products for consumers. In such states or countries some exonerations or limitation of this warranty may not be applied to the user. This warranty is subject to modifications without any prior notice.

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