ACTIVESHOP®

USER MANUAL

Sillon Luxury electric beauty chair, 3 motors, foot control

C502 3 motors / 150223



Importer:
ACTIVESHOP SP. Z O.O.
ul. Graniczna 8b, bud. DC2A
54-610 Wrocław, Poland
www.activeshop.com.pl
Rev. A/02/2026

Scan for different languages

Thank you for choosing our product.
Please read the user manual before first use and keep it for future reference.

The **Sillon Luxury beauty chair** combines modern design, intuitive operation, and wide functionality. Equipped with three motors controlled by a wired remote control, it allows for precise adjustment to individual preferences. The foot drive provides the user with comfort and precision in adjusting the position of the chair without using their hands.

Technical data	
Input voltage	220-230 V
Electrical power	200W
Rated frequency	50Hz
Upholstery material	Synthetic leather
Control	Wired remote control and foot control
Room temperature	(+)10°C – (+)40˚C
Ambient humidity	30% - 75

Note


Cardboard and outer boxes are used for packaging, with a total of 10 layers of paper to ensure proper protection of the product. We recommend keeping the original packaging for future use, especially in case of returns during the warranty period.

The chair can only be lifted from the headrest and footrest positions, not from both sides of the armrests. For safety reasons, two people are needed to remove the chair from the cardboard box.

The electric beauty chair must be used in closed, dry rooms. Temperature
The room temperature should be between (+)10 °C and (+)40 °C, and the relative humidity between 30% and 75%. If this beauty chair has been exposed to temperatures below 0 °C, e.g. during transport, wait at least 5 hours before starting it up to allow the device to reach room temperature.

Device construction:

1. Headrest
2. Backrest
3. Armrest
4. Seat
5. Footrest
6. Wired remote control
7. Foot control
8. Lifting cylinder
9. Plastic cover

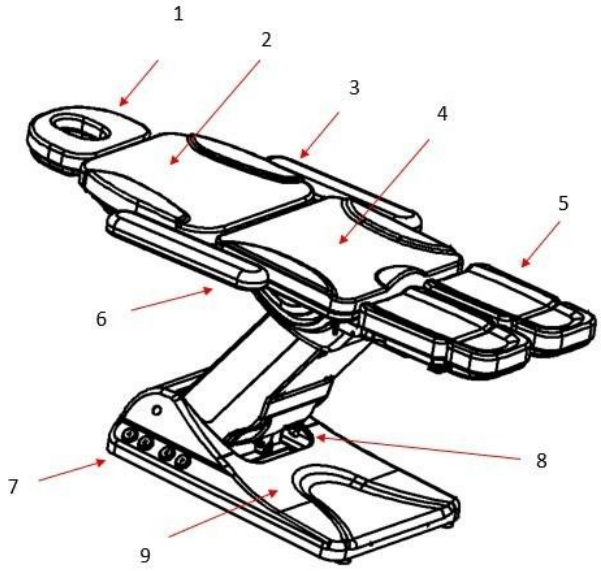


Chair assembly instructions

- Place the chair in the location where it will be used
- Plug the power cord into a power outlet


Description of remote control functions

1. **Seat height adjustment.**

   The height adjustment range is 60 cm – 82 cm from the floor

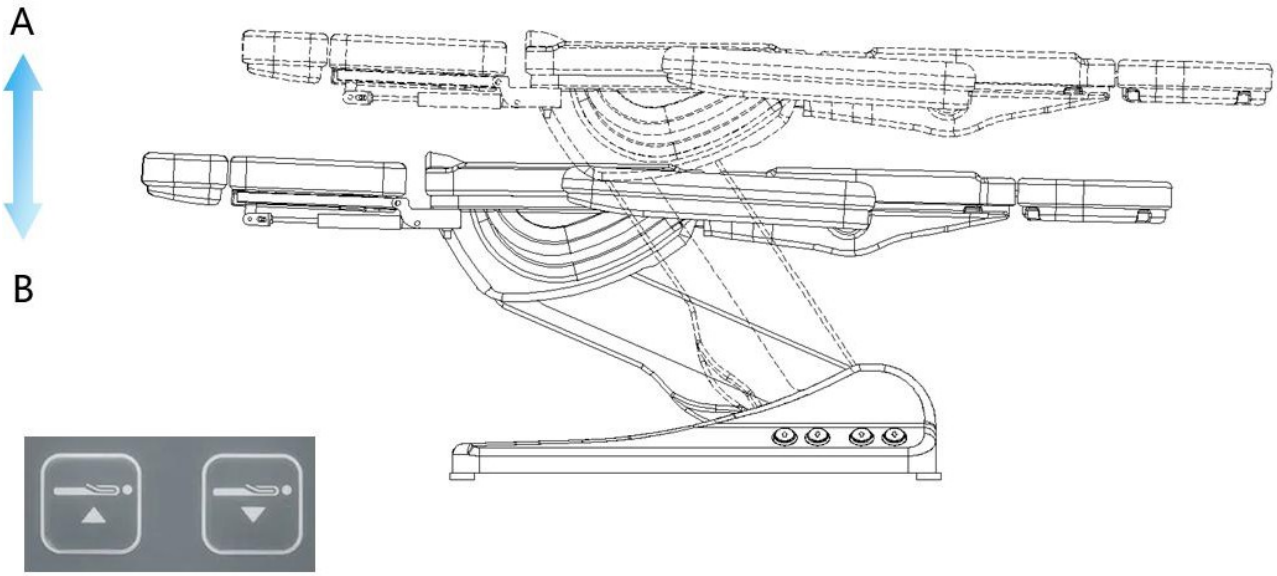


2. **Backrest angle adjustment –** range 0–80°

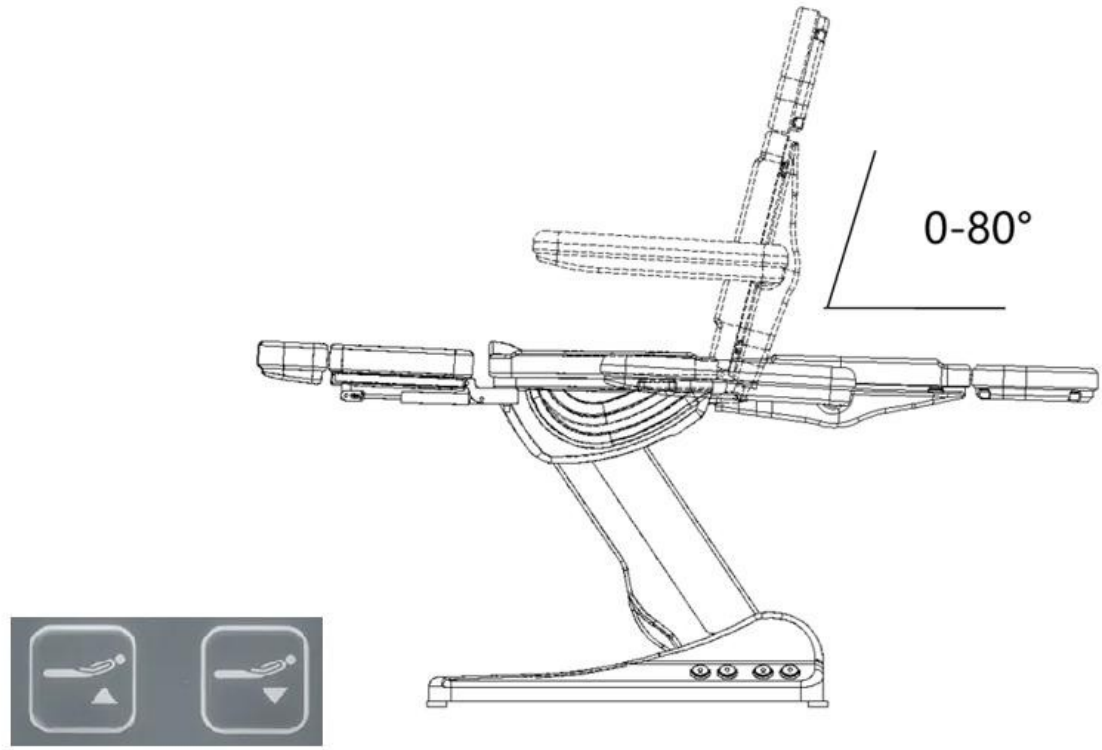

3. **Leg angle adjustment –** range 0–90°

   Press the pneumatic rod handle to change the footrest position.

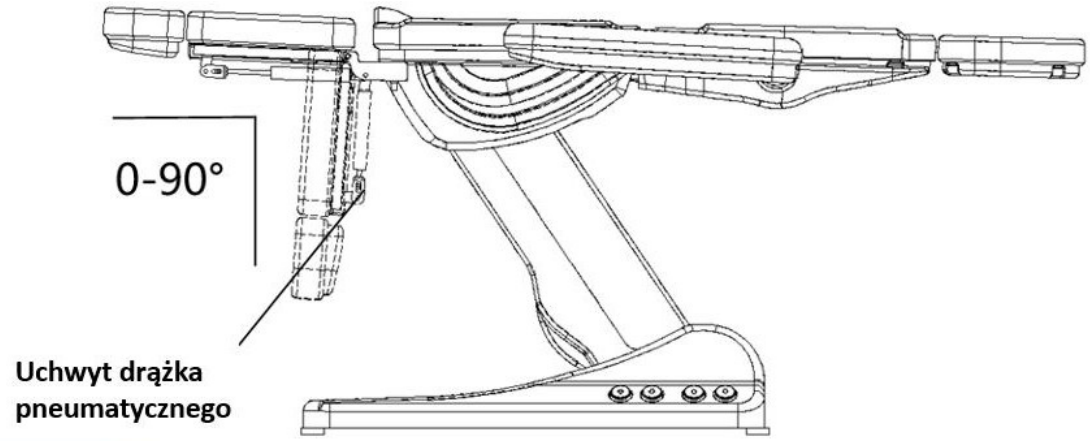




4. **Base bottom lighting**



5. **Note: The heating function is not active in this model. The remote control is universal and can be used with many models.**

Telescopic headrest and footrest extension function – 12 cm extension

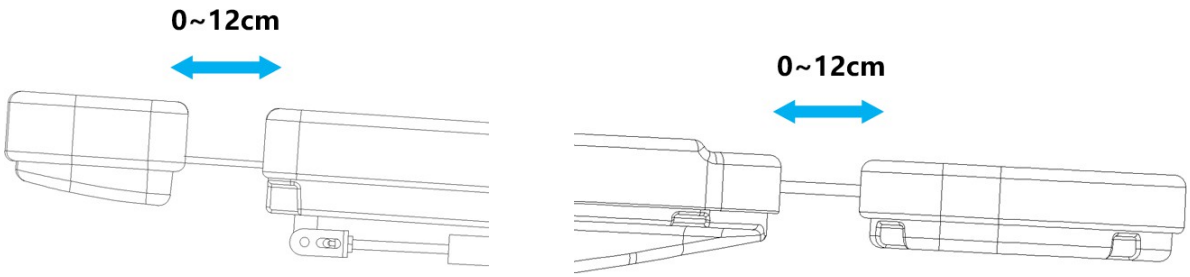

Headrest


Footrest

Armrest installation

The chair must be in a horizontal position. Insert the armrest into the hole on the side of the backrest, move it to the right and left until it clicks into place. Do the same with the other armrest.
The armrest is perpendicular to the ground regardless of the angle of the seat. Its range of motion is 180°.

Note

The maximum load is 150 kg. Please do not sit directly on the reclined backrest, armrests
and footrest.

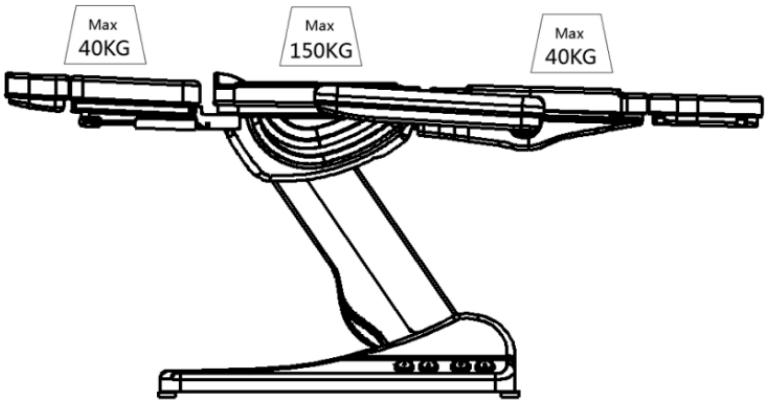

Maintenance

To keep the surface of the beauty bed in good condition, it should be cleaned and maintained regularly.

Before cleaning, remember to unplug the power cord from the power outlet.

- Metal and plastic surfaces and controls should be cleaned with a damp cloth and a mild cleaning solution. Use a small brush in corners and other hard-to-reach areas. We recommend not using too much liquid when cleaning.
- Stainless steel surfaces - clean with a mild detergent solution. Ammonia and most solvents can be used to remove difficult and stubborn stains. Avoid chlorine-based solutions.
- Painted or chrome-plated metal surfaces - do not use solvents or acids for cleaning.
- Bed surface - use a clean, soft cloth. A small amount of a special cleaning and impregnation agent for synthetic leather can be used.

All surfaces must be dry before connecting to the power supply and reuse. Do not clean at

high temperatures and high humidity.


Do not expose the device to factors that may lead to contamination with liquids.

Recommendations:

1. Ensure that you have installed the device on a flat and stable surface to reduce the risk of falling and damage.
2. Never leave the device unattended while it is plugged in.
3. Use the device in accordance with the intended use described in the manual.
4. Never use this device if the cord or plug is damaged.
5. Do not use the device if it is not working properly or if it has been dropped, damaged, or fallen into water.
6. Do not pull directly on the cord to remove the plug from the outlet.
7. This product is for indoor use only. Do not expose it to direct weather conditions.

8. Clean the device periodically to prevent dirt accumulation. It is best to clean the device with a damp cloth.
9. Do not wrap the power cord around the device.

Device dimensions

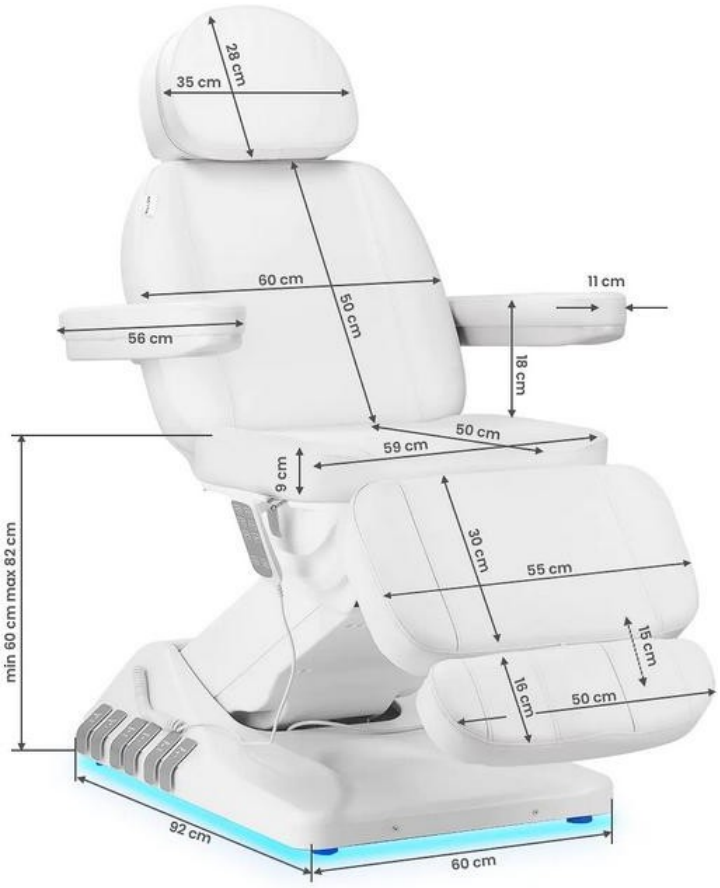

Safety precautions:

1. The device is professional equipment for use in beauty salons. Only qualified beauty personnel should operate the device to avoid undesirable effects during treatment.
2. Do not disassemble the device or attempt to use it for purposes other than those specified in this user manual. All service and repair work must be performed exclusively by an authorized service technician (ACTIVESHOP Sp. z o.o.).
3. Do not install or operate the device near liquids, wet areas, or with wet hands. Ensure that no liquids get into the device.
4. Immediately unplug the device and contact an authorized service center provided by ACTIVESHOP Sp. z o.o. in the following cases: a) liquid has entered the device,
   b) the device emits a strange odor, smoke, or strange noises, c) the power cord has been damaged, d) the device has been knocked over or its casing has been damaged.
5. Do not place any objects on the cables. Do not install the device in places where someone could step on the power cord.
6. For safety reasons, unplug the power cord from the outlet after each use of the device.
7. Do not place foreign objects inside the device. If an unwanted object gets inside, unplug the power cord from the power outlet and contact your distributor or manufacturer.
8. Do not place the device on an unstable cart, shelf, or vehicle. If the device falls,

may be damaged.

9. In order to guarantee the specified service life of the device, it is recommended to wait one minute before pressing the power button again (between switching on and off and vice versa).
10. The power supply must comply with the specifications on the rating plate. Otherwise, there is a risk of incorrect operation, damage, or even burning of the device.

CAUTION

1. In the event of damage to the power cord or other malfunctions, contact an authorized service center to avoid danger.
2. This device can be used by children aged 8 years and older, by persons with reduced physical and mental abilities, as well as by persons unfamiliar with the operation of this equipment, if they are provided with appropriate supervision or instruction on how to use it safely and understand the risks involved.
3. Do not leave the device unattended when children or persons who are not familiar with its use are in the vicinity.
4. The device is not intended for use by children. Do not allow children to play with the device.
5. This electrical device is not a toy. Incorrect use or contact with an electrical device may pose a risk to health or lead to death, so the device must be used and stored out of the reach of children.

Disposal of the device



In order to protect the environment, in accordance with Directive 2012/19/EC on waste electrical and electronic equipment, used electrical equipment must be collected separately and taken to collection points for waste electrical and electronic equipment. Follow the rules for disposing of this type of equipment and do not dispose of it in household waste containers.

At the end of their useful life, these devices still contain valuable recyclable components, so they should be sent for recycling to appropriate collection points or to ACTIVESHOP Sp. z o.o.

Storage and transport

The device should be transported at a temperature between (–)20°C and (+)40°C only in its original packaging. It should be stored in a well-ventilated room with a relative humidity of no more than 90%. Operating temperature from (+)10°C to (+)40°C. During processes storage and transport, do not throw or drop the device, expose it to mechanical shocks or atmospheric factors.