moretti®

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EN: USED ONLY WITH DETACHABLE CONTROLLER SSO3. A TOO LONG APPLICATION CAN RESULT IN SKIN BURNS. DO NOT DRY CLEAN. DO NOT USE WITH INFANTS, A HELPLESS PERSON OR A PERSON INSENSITIVE TO HEAT. HAND WASH ONLY. DO NOT USE IF WET.

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EN

1. CODES

LTK700 Multi-purpose heating pad - 9 temperatures
LTK710 Neck and back heating pad - 9 temperatures
LTK720 Lumbar heating belt - 3 temperatures
LTK730 Neck and shoulders heating cape - 3 temperatures
LTK740 Neck and back heating cape - 3 temperatures
LTK750 Back heating cape - 3 temperatures
LTK770 Electric feet warmer - 3 temperatures

2. INTRODUCTION

Thank you for choosing a Moretti heating pad.

Before using the product, please read this manual carefully, especially the safety instructions, and follow them during use. Keep the manual for future reference and pass it on together with the product if it is transferred to a third party.

Inside, you will find simple instructions for proper use and useful advice to ensure your safety.

For any questions or assistance, please contact your trusted retailer.

Note: Check that all parts of the product have not been damaged during the shipment. In case of damage do not use the product and contact the retailer for further instructions.

3. REGULATORY COMPLIANCE

This product complies with the essential requirements set out by the applicable European legislation for low-voltage electrical devices.

In particular:

- Compliant with the European Low Voltage Directive (LVD) 2014/35/EU.
- Product classified as a non-medical device.
- Designed and tested in accordance with the following technical standards:
 - **EN 60335-1:2012** (Safety of household electrical appliances)
 - **EN 60335-2-17:2013** (Particular requirements for heating pads and electric blankets)
 - **EN 62233:2008** (Electromagnetic fields)

The maximum temperature delivered by the device is 65°C ± 5°C, in compliance with the limits set by the relevant standards.

4. SAFETY WARNINGS

4.1 Electrical hazards

- Before connecting the appliance, check that the voltage indicated on the rating label matches the voltage of the electrical supply.
- DO NOT operate the controller with wet hands.
- DO NOT pull the power cable to unplug the appliance.
- DO NOT immerse the power cable or the remote controller in water.
- Always unplug the appliance before cleaning, performing any maintenance, or when the device is not in use.
- In the event of a fault or malfunction, switch off the device and do not attempt to repair it.
- Contact an authorised service centre.
- If the cable is damaged, it must be replaced by authorised personnel.
- Use only the controller indicated on the product label.

4.2 Use with children and vulnerable persons

- NOT suitable for children under 3 years of age.
- Children must always be supervised.
- Do not use on infants, sick individuals, or persons insensitive to heat.
- Persons with reduced physical, sensory, or mental capabilities should use the device only under supervision.
- Do not use on individuals with pacemakers.
- This product is not intended for medical or hospital use.

4.3 General safety and packaging

- After opening, check that the product is not damaged.
- In case of damage, do not use it and contact your retailer.
- Keep packaging materials out of reach of children: risk of suffocation.

4.4 Main prohibitions and improper use

- Use only as a heating pad for domestic use.
- DO NOT use the device if wet or in humid environments (e.g. bathroom).
- DO NOT use while folded.
- DO NOT place the controller or power supply under or on top of the product during use.
- DO NOT leave objects or clothing inside when the product is folded.
- Allow to cool before storing.
- Regularly check the condition of the product; do not use it if worn or damaged.
- DO NOT use to warm animals.
- DO NOT use in the presence of flammable gases or near gas sources.
- DO NOT place combustible materials on the product (paper, plastic, cardboard, etc.).
- DO NOT expose to heat sources, frost, direct sunlight, or moisture.
- DO NOT insert needles or sharp objects.

⚠ **WARNING!**
There is a risk of injury in case of improper use.

5. FIRST USE

1. Connect the controller.
 Insert the controller connector into the heating pad socket.
2. Plug into the power supply.
 Insert the plug into an electrical outlet.
3. Turn on the device.
 Select the temperature level:
 - 1 = low level, for continuous use
 - 2–5 = medium level
 - 6–9 = high level
4. Position the heating pad
 Place it on the area of the body to be treated (back, neck, abdomen, etc.).

5.1 Operation

- The device automatically switches off after 90 minutes.
- To restart: set to 0, wait a few seconds, then select the temperature again.

5.2 Switching off

1. Set the selector to “0”.
2. Unplug from the electrical outlet.

EN

WARNING (IMPORTANT)
- Do not use if wet.
- Do not fall asleep with the device switched on.
- Do not use on children or persons insensitive to heat.
- Do not fold during use.

5.3 Cleaning (essential)
- Always disconnect from the power supply.
- Hand wash only (max 30°C).
- Do not use bleach.
- Allow to dry completely before use.

6. SYMBOLS

Product code

CE mark

Month and year of production

Instruction manual

Disposal

Insulation class II

Do not use bleach

Hand wash


Do not dry-clean

Do not dry the heating pad in the dryer

Do not iron


7. GENERAL DESCRIPTION

7.1 LTK700-LTK710

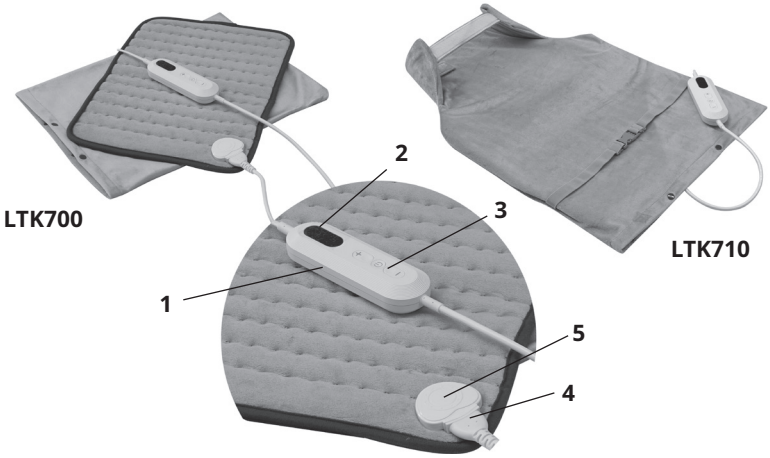Main components
1. Remote control
2. Temperature selector
3. Display
4. Connector
5. Connector base

7.2 LTK720

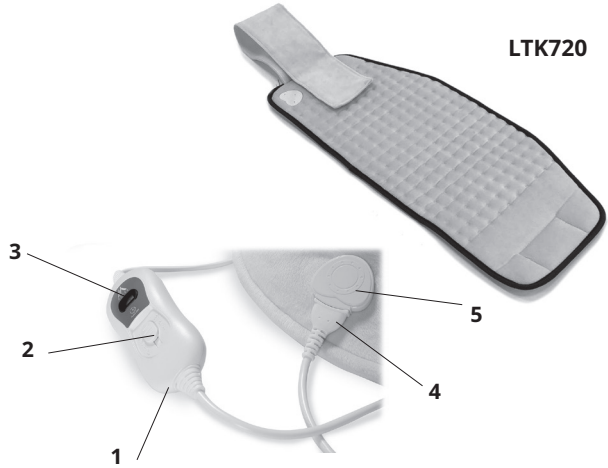
Main components
1. Remote control
2. Temperature selector
3. Display
4. Connector
5. Connector base

Image for illustration only


EN

7.3 LTK730 - LTK740 - LTK750

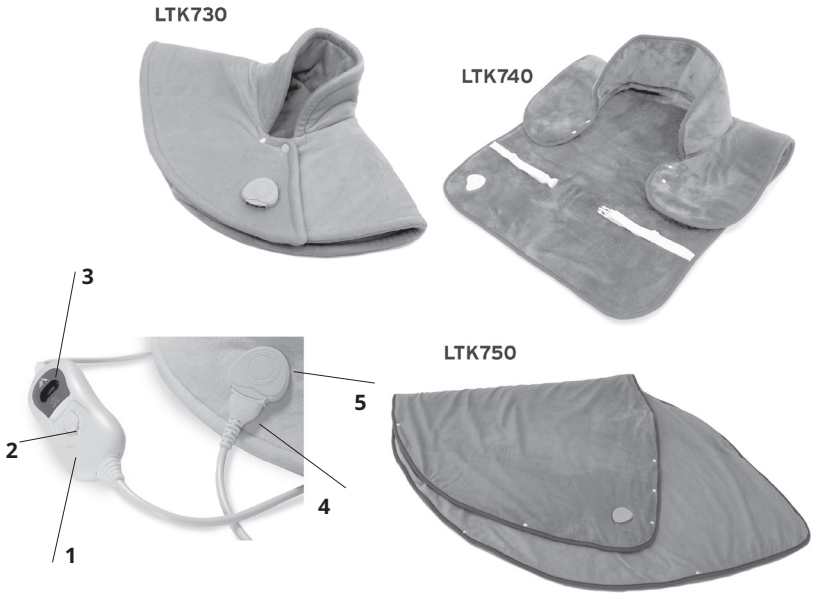Main components
1. Remote control
2. Temperature selector
3. Display
4. Connector
5. Connector base

7.4 LTK770

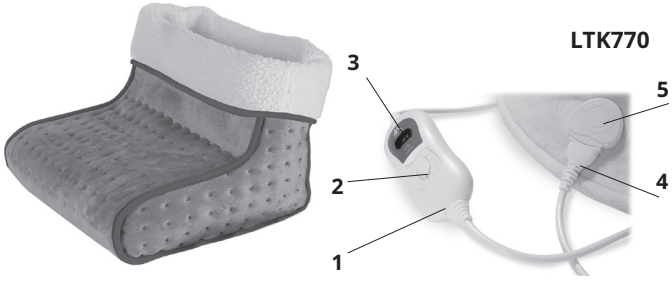Main components
1. Remote control
2. Temperature selector
3. Display
4. Connector
5. Connector base

Image for illustration only

8. INSTRUCTIONS FOR USE

The heating pad is designed to provide localized heat to the body. It can be used on various areas, including: neck (cervical area), shoulders, back, lower back, abdomen, joints, hands, and feet, depending on the model.

WARNING!
- Use only with the controller supplied with the product.
- The controller model is indicated on the product label.

9. USE

1. Connection

Insert the controller connector (4) into the socket (5) of the heating pad.

2. Power on

Plug the device into the electrical outlet.

3. Temperature selection depending on the model available.

Turn the selector to set the desired level:

Level	Function
0	Off
1	Low level, for continuous use
2–5	Medium heat
6–9*	High heat

* available only on 9-level models

10. CLEANING AND MAINTENANCE

WARNING! BEFORE STARTING
- Disconnect the heating pad from the power supply.
- Unplug the connection cable.
- Ensure the product is completely cool.

10.1 Washing

Hand wash only; Maximum temperature: 30°C

Do not use bleach

Do not dry-clean

- Use a mild detergent.
- Do not use chemicals or abrasive products.

10.2 Drying

Do not dry in the dryer

- Allow to air dry.
- Do not place on radiators or heat sources.
- Do not speed up drying with electrical devices.

Image for illustration only

Do not iron

EN

10.3 Before reuse
Ensure the heating pad is completely dry before use.

11. GENERAL CONDITIONS OF DISPOSAL

The product is packaged in recyclable materials. Dispose of it in compliance with the environmental protection standards.

The product must be disposed of appropriately, in compliance with Directive 2012/19/EU (WEEE). The recyclable materials contained in the appliance must be recovered so as to prevent environmental pollution. For further information, contact the local disposal body or the appliance dealer

12. TECHNICAL FEATURES

moretti
CON TE, CON CURA.

NO BLEACH NO IRONING NO DRY CLEAN NO WASHER DRY

Model:
Voltage: **220-240 ~ 50Hz**
Rated power input: **100W**

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13. TROUBLESHOOTING

WARNING!
If the electrical device does not work properly and it is defective, do not attempt to repair it by yourself.

PROBLEM	POSSIBLE CAUSE	SOLUTION
The device does not work.	The device is not receiving power.	Check if: **a.** Check that the plug is properly inserted; **b.** Verify the cable connection; **c.** Check the electrical outlet; **d.** Check the fuse.
	No temperature level has been selected. The temperature selector is in position 0.	Move the temperature selector
	The protection against overheating has been activated.	Have the device checked by the manufacturer.

14. WARRANTY

Every Moretti product is guaranteed for 2 (two) YEARS from the purchase date, except the restrictions specified below. The warranty shall not apply in case of damage caused by improper use, abuse, alterations of the product or any other use not in compliance with the user manual. The proper use of the product is indicated in the user manual. Moretti S.p.A. shall not be liable for any damages or personal injuries due to incorrect assembly/improper use of the device, not compliant with the user manual. Moretti does not guarantee its products against damages or faults in the following cases: natural disasters, non-authorized maintenance or reparations, damages caused by electrical supply problems (if provided), use of parts or components not provided by Moretti, failure in accomplishing the guidelines and instructions for use, non-authorized alterations, damages due to the transport (different from Moretti original delivery) or failure in maintaining the device as indicated in the manual. The components subject to wear are not covered by this warranty if the damage is due to the normal use of the product.

15. EXEMPT CLAUSES

Except what clearly specified in this warranty and in accordance with the law, Moretti S.p.A. does not offer its customers any other declaration, warranty or condition, expressed or implicit, including any declaration, warranty or condition of merchantability, of non violation or non interference, of suitability for a specific aim. Moretti S.p.A. does not guarantee that the use of Moretti products will be continuous and faultless. The duration of possible implicit warranties that may be applied in accordance with the law is limited to the warranty period. Some states or countries do not allow limitations to the duration of an implicit warranty, or the exoneration, or the limitations of accidental or indirect damages related to products for consumers. In such states or countries some exonerations or limitation of this warranty may not be applied to the user. This warranty is subject to modifications without any prior notice.

Please consult our website for the most up-to-date version of the user manual.

Moretti S.p.A.
Via Bruxelles 3 - 52022 Cavriglia (AR) - Italy
www.morettispa.com - info@morettispa.com
Made in P.R.C.

Moretti S.p.A.
Via Bruxelles 3
52022 Cavriglia (AR) Italy

info@morettispa.com
www.morettispa.com